

JRC TECHNICAL REPORTS

Pyrrolizidine alkaloids in herbal tea and honey

*Report on the 2017
Proficiency testing
scheme*

Breidbach, A.
Tamošiūnas, V.

2017



This publication is a Technical report by the Joint Research Centre (JRC), the European Commission's science and knowledge service. It aims to provide evidence-based scientific support to the European policymaking process. The scientific output expressed does not imply a policy position of the European Commission. Neither the European Commission nor any person acting on behalf of the Commission is responsible for the use that might be made of this publication.

Contact information

Name: Andreas Breidbach
Address: Retieseweg 111, 2440 Geel, Belgium
Email: andreas.breidbach@ec.europa.eu
Tel.: +32 (0)14 571205

JRC Science Hub

<https://ec.europa.eu/jrc>

JRC109959

EUR 28960 EN

PDF ISBN 978-92-79-77167-5 ISSN 1831-9424 doi: 10.2760/219405

Luxembourg: Publications Office of the European Union, 2017
© European Union, 2017

Reuse is authorised provided the source is acknowledged. The reuse policy of European Commission documents is regulated by Decision 2011/833/EU (OJ L 330, 14.12.2011, p. 39).

For any use or reproduction of photos or other material that is not under the EU copyright, permission must be sought directly from the copyright holders.

How to cite this report: Breidbach A, Tamošiūnas V, *Pyrrrolizidine alkaloids in herbal tea and honey—Report on the 2017 Proficiency testing scheme*, EUR 28960 EN, Publications Office of the European Union, Luxembourg, 2017, ISBN 978-92-79-77167-5, doi 10.2760/219405, JRC109959.

All images © European Union 2017, except: Title By I, Tony Wills, CC BY 2.5,
<https://commons.wikimedia.org/w/index.php?curid=3323872>



European Union Reference Laboratory
Mycotoxins



268-PT Accredited by the
Belgian Accreditation Body (BELAC)

Executive Summary

Pyrrrolizidine alkaloids (PAs) and their N-oxides (PANOs) are plant toxins which can enter the food chain through different paths. Two affected foods are herbal infusions and honey. This proficiency testing scheme was executed to assess the capabilities of laboratories to determine PAs. 29 laboratories from nine EU Member States plus Singapore registered. On 04. and 06.09.2017 test items and documentation were dispatched to all of those laboratories.

By the dead line of 24.10.2017 26 laboratories had reported back results and filled in a questionnaire. Test item HO (acacia honey) was fortified with six PAs/PANOs (Echimidine, Integerrimine, Intermedine, Senecionine, Seneciphylline-NO, and Senkirkine) and 23 laboratories reported results for this item. The same number of laboratories reported for test item HT (herbal infusion) which was naturally contaminated with four PAs after extraction under reductive conditions (Integerrimine, Retrorsine, Senecionine, and Senecivernine). Laboratories had to report the sums of PA and its respective PANO.

Satisfying outcomes could only be registered for Senecionine in test item HT and for Echimidine, Intermedine, and Senkirkine in test item HO with 74 %, 85 %, 85 %, and 91 %, respectively, of reported results having a z'-score smaller or equal to $|2|$. Only four laboratories reported for Integerrimine in both test items. Contrary to test item HT, Senecionine analysis in test item HO showed very unsatisfactory results. Of the 22 z'-scores calculated for Senecionine nine (41 %) were larger than 3. Senecivernine measurements in test item HT showed a similarly unsatisfying outcome with 47 % of reported results having z'-scores larger than 3.

Only three laboratories out of the 26 were able to test for all 10 measurands and only one reported all 10 values with z'-scores smaller or equal to $|2|$. Overall only five laboratories obtained satisfactory z'-scores ($\leq |2|$) for all their reported results. There are two groups of three isomeric PAs/PANOs each which apparently caused, for a number of laboratories, problems with quantification. This is an issue which deserves heightened attention.

The questionnaire contained queries regarding accreditation and experience, preparation conditions for the two test items, chromatographic separation conditions, detection conditions, calibration approach, and a comments section. The answers were evaluated and for selected questions their correlation to the z'-score of Senecionine in test item HO or Senecivernine in test item HT was analysed. For none of the tested questions a significant influence could be shown.

Table of Contents

| | |
|---|----|
| Executive Summary..... | 2 |
| 1. Introduction | 4 |
| 2. Scope..... | 4 |
| 3. Set-up of PT scheme | 5 |
| 3.1. Time frame | 5 |
| 3.2. Confidentiality..... | 5 |
| 3.3. Distribution | 5 |
| 3.4. Instructions to participants..... | 5 |
| 3.5. Standard deviation for proficiency assessment..... | 5 |
| 4. Test items..... | 6 |
| 4.1. Preparation | 6 |
| 4.2. Homogeneity..... | 6 |
| 4.3. Stability | 6 |
| 5. Assigned values and their uncertainties | 6 |
| 6. Evaluation of the results | 7 |
| 6.1. Test item HT | 7 |
| 6.2. Test item HO | 8 |
| 6.3. "Smaller than" / "Larger than" values and uncertainties | 10 |
| 6.4. The questionnaire | 11 |
| 7. Conclusions | 13 |
| 8. Acknowledgements..... | 13 |
| References: | 13 |
| Annex 1: Instructions for PT scheme | 14 |
| Annex 2: Proficiency testing materials receipt form | 17 |
| Annex 3: Homogeneity testing | 19 |
| Annex 4: Stability testing | 20 |
| Annex 5: Determination of the assigned values and their uncertainties | 22 |
| Annex 6: Graphs..... | 26 |
| Annex 7: Questionnaire data | 32 |
| Annex 8: Participating laboratories | 44 |

1. Introduction

Pyrrolizidine alkaloids (PAs) are secondary metabolites of plants used as defence against plant eating animals. It is estimated that over 6000 plant species worldwide, mainly from the botanical families Boraginaceae (e.g. *Heliotropium* spp.), Asteraceae (e.g. *Senecio* spp.) and Fabaceae (e.g. *Crotalaria* spp.), produce at least 600 different PAs. Those with an unsaturated bond in position 1,2 of the pyrrolizidine ring system (Figure 1) are hepatotoxic. Additional diversity is created due to oxidation of the pyrrolizidine nitrogen to the N-oxide (PANO) which doubles the possible number of PAs.

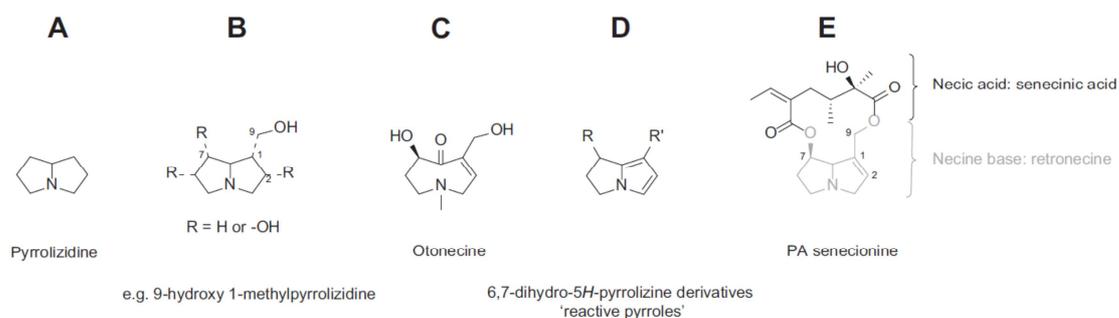


Figure 1: Structural features of PAs. (A) core structural motif pyrrolizidine (1,2,3,6,7,8-hexahydro-5H-pyrrolizine); (B) general description of the main necine base parts of naturally occurring PAs including the common necine base numbering; (C) necine base otonecine; a core structural motif of otonecine-type PAs; (D) general pyrrolizidine structure motif and (E) structural example of 1,2-unsaturated ester PA senecionine (Figure taken from [1])

PAs may enter the food chain via different routes. Two possibly affected end products are teas / herbal infusions, if PA containing weeds had contaminated the leaves / herbs, and honey, if PA containing pollen was collected by the bees. In 2016 the European Food Safety Authority (EFSA) published a new scientific report about dietary exposure to PAs [1]. In that report teas / herbal infusions were identified as the largest contributors to total PA exposure. It has been recommended that for relevant food commodities, as honey and teas / herbal infusions, the effort to collect PAs occurrence data should be continued.

In this PT scheme a naturally contaminated rooibos herbal blend for infusion and a spiked acacia honey were sent out to 29 laboratories for PAs analysis. The two test materials were prepared, tested for homogeneity, stability, and characterized with respect to contamination level in-house. This report presents the results of these tests and characterizations as well as the reported results from the participating laboratories.

2. Scope

This PT scheme was executed to assess the capabilities of the participating laboratories to determine one or more of the following PAs and/or their corresponding PANOs:

Echimidine, Erucifoline, Europine, Indicine, Integerrimine, Intermedine, Jacobine, Lasiocarpine, Lycopsamine, Monocrotaline, Retrorsine, Senecionine, Seneciphylline, Senecivernine, Senkirkine, Trichodesmine,
in an herbal blend for infusion and acacia honey.

3. Set-up of PT scheme

3.1. Time frame

The PT registration was open from 19.05.2017 to 15.08.2017. Registration could be done through a dedicated web site. The test materials were dispatched to the participants on 04.09. and 06.09.2017. All shipments reached their respective recipients within 2 days. The initial deadline for reporting the results was 20.10.2017. This deadline was eventually extended to 24.10.2017 because two laboratories reported problems with entering the results.

3.2. Confidentiality

The procedures used for the organisation of PTs are accredited according to ISO 17043:2010 [2] and guarantee that the identity of the participants and the information provided by them is treated as confidential.

3.3. Distribution

In all, 29 laboratories from nine EU Member States plus Singapore registered for participation. The following was sent to each of them:

- one container with 40 – 50 g of a herbal blend for infusion naturally contaminated with PAs;
- one container with 50 – 70 g of acacia honey spiked with PAs
- the "Instructions for PT scheme" letter with the individual reporting password (Annex 1: Instructions for PT scheme);
- the "Proficiency testing materials receipt form" (Annex 2: Proficiency testing materials receipt form);
- a "Submission of results" document explaining the use of the reporting interface.

3.4. Instructions to participants

The "Instructions for PT scheme" letter included with the test item shipment gave detailed instructions to the participants. In brief, the participants were informed that one or more of the following PAs:

Echimidine, Erucifoline, Europine, Indicine, Integerrimine, Intermedine, Jacobine, Lasiocarpine, Lycopsamine, Monocrotaline, Retrorsine, Senecionine, Seneciphylline, Senecivernine, Senkirkine, Trichodesmine,

could be present in either one of the two test items. The participants had to check the test items for integrity and store them upon receipt at temperatures between -15 °C and -25 °C.

Instructions on how to access the reporting web site and how to enter results were also included.

3.5. Standard deviation for proficiency assessment

Scoring of the reported results was to be done acc. to ISO 13528:2015 and was to be based on an independently determined reference value (x_{PT}) for each PA and a target standard deviation σ_{PT} of 0.22 x_{PT} acc. to the following equation:

$$z = \frac{(x_i - x_{PT})}{\sigma_{PT}}$$

with x_i the reported result of an individual PA from a laboratory. No scoring was to be provided for reported sums of total PAs.

For those analyte contents for which the associated uncertainty of x_{PT} ($u(x_{PT})$) was larger than $0.3 \sigma_{PT}$ the following equation was used:

$$z' = \frac{(x_i - x_{PT})}{\sigma'_{PT}}$$

with the adjusted target standard deviation $\sigma'_{PT} = \sqrt{\sigma_{PT}^2 + u^2(x_{PT})}$.

4. Test items

4.1. Preparation

The herbal infusion blend (test material HT, mostly rooibos) was purchased from a provider within the EU. It was milled with a ZM200 centrifugal mill with a 0.5 mm sieve (Retsch, Hahn, Germany). After thorough homogenization the material was filled into 96 plastic screw-capped containers.

The acacia honey (test material HO) was courtesy of Breitsamer + Ulrich GmbH & Co KG (München, Germany). The honey was practically free of PAs (none detected above 3 $\mu\text{g}/\text{kg}$) and spiked with Echimidine, Integerrimine, Intermedine, Senecionine, Seneciphylline-NO, and Senkirkine. It was then thoroughly mixed with a ploughshare mixer (Lödige, Paderborn, Germany) and filled into 73 plastic screw-capped containers.

4.2. Homogeneity

For testing the homogeneity ten containers of HT and eight containers of HO were selected. Test portions were extracted/cleaned-up with SPE acc. to an internal protocol and measured in randomized order with a LC-HRMS. The obtained signals, without conversion to concentration units, were evaluated for sufficient homogeneity acc. to ISO 13528 Annex B [3]. More details can be found in Annex 3: Homogeneity testing.

For all analytes in both test materials the results indicated sufficient homogeneity.

4.3. Stability

The stability of the analytes in the two test materials was determined through an isochronous study. For this purpose randomly selected test units were stored frozen (-15 to -25 °C) and at 4 °C, and 40 °C for up to eight weeks. At the day of measurement all test units were prepared and measured as described for the homogeneity testing. More details can be found in Annex 4: Stability testing.

All analytes in both test materials were sufficiently stable during the PT period.

5. Assigned values and their uncertainties

For the scoring of this PT scheme individual mass fractions for all the PAs in the two test materials were determined. This was done with gravimetric standard addition using an internal standard

(ISTD) acc. to Hauswaldt et al [4]. These measurements were executed in early October for test item HT and early November for test item HO. This was due to time constraints and happened without prior knowledge of reported results. More details can be found in Annex 5: Determination of the assigned values and their uncertainties.

The assigned values x_{PT} and their associated standard uncertainties $u(x_{PT})$ were determined acc. to ISO 13528 clause 7 [3] and are listed in Table 1 for test item HT and Table 2 for test item HO.

6. Evaluation of the results

6.1. Test item HT

Test item HT was naturally contaminated with four PAs/PANOs at detectable levels: Integerrimine, Retrorsine, Senecionine, and Senecivernine. Table 1 lists relevant parameters. The reference values (x_{PT}) are for the sum of PA/PANO after reductive extraction. Since the associated uncertainties of the reference values were larger than $0.3 \sigma_{PT}$ in all cases the adjusted target standard deviation σ'_{PT} was used for scoring. The robust mean (\bar{x}_{rob}) was computed from the reported results using Algorithm A [3, Annex C] following the Recommendation 1 in [5]. A robust mean could not be calculated for the reported results of Integerrimine because of too few results and for Senecivernine because of too high a contribution of minor modes to the multi-modal distribution of reported results. For the other two robust means, the one of Senecionine falls within the expanded uncertainty range ($k=2$, ~ 95 % confidence) around x_{PT} , the one of Retrorsine does not.

For the four PAs present in test item HT a satisfying outcome can only be attested for the analysis of Senecionine. It was present at the highest level and 17 out of 23 (74 %) reported results showed z' -scores $\leq |2|$. No reported result had a z' -score $> |3|$. At the opposite end is the outcome for the measurements of Senecivernine with 47 % of reported results with z' -scores larger than 3. For Integerrimine only four laboratories reported a result, two of them as "smaller than", even though it was present at the second highest level in this naturally contaminated material. Table 3 is listing all reported results and the corresponding z' -scores. The figures in Annex 6: Graphs for Test item HT visualise these results.

Integerrimine, Senecionine, and Senecivernine are three isomeric PAs with the elemental formula ($C_{18}H_{25}NO_5$). This makes these three compounds practically indistinguishable for mass spectrometry which was the detection mode used by all participating laboratories. Only proper chromatographic separation, which is challenging, can provide reliable results for the individual compounds. Insufficient separation between Integerrimine and Senecivernine, and false identification may have led to the large number of too high results for Senecivernine. A compounding factor may be that the majority of participating laboratories were apparently not testing for Integerrimine and, therefore, might not have been aware of any coelution issues. Integerrimine is not included in the list of PAs provisionally selected by the European Commission [1].

Laboratory 14 reported to have also detected Senkirkine at 10 µg/kg. No other laboratory reported any other PA above its respective LOQ. Laboratory 3 reported the sum of all PAs as "retronecine equivalents" at 78 µg/kg.

Table 1 Parameters of the four PAs and summary of the outcome in test item HT

| | Integerrimine | Retrorsine | Senecionine | Senecivernine |
|-------------------------------|---------------|------------|-------------|---------------|
| x_{PT} [µg/kg] | 30.6 | 19.8 | 80.4 | 5.7 |
| $u(x_{PT})$ [µg/kg] | 2.8 | 1.5 | 7.8 | 1.2 |
| \bar{x}_{rob} [µg/kg] | | 15.4 | 65.4 | |
| σ_{PT} [µg/kg] | 6.7 | 4.3 | 17.7 | 1.2 |
| $u(x_{PT}) > 0.3 \sigma_{PT}$ | YES | YES | YES | YES |
| σ'_{PT} [µg/kg] | 7.3 | 4.6 | 19.3 | 1.7 |
| Number reported | 4 | 20 | 23 | 15 |
| "<" or ">" | 2 | 7 | 1 | 2 |
| $z' \leq 2 $ | 2 | 9 | 17 | 6 |
| $ 2 < z' \leq 3 $ | 0 | 2 | 5 | 0 |
| $z' > 3 $ | 0 | 2 | 0 | 7 |

6.2. Test item HO

Test item HO was an acacia honey free of PAs. It was spiked with Echimidine, Integerrimine, Intermedine, Senecionine, Seneciophylline-NO, and Senkirkine. Table 2 lists relevant parameters. These parameters were determined as described for test item HT above. Again, for Integerrimine too few results were reported to calculate a robust mean. Of the other robust means, the ones for the content of Echimidine and Senkirkine lay within the expanded uncertainty ranges ($k=2$, ~ 95 % confidence) around the respective x_{PT} , the ones for Intermedine, Senecionine, and Seneciophylline were below the lower limit of the expanded uncertainty ranges. Also, as for test item HT, the associated uncertainties of the reference values were larger than $0.3 \sigma_{PT}$ in all cases but for Senecionine. Yet for reasons of consistency the adjusted target standard deviation σ'_{PT} was used for scoring of all PAs including Senecionine.

Table 2 Parameters of the six PAs and summary of the outcome in test item HO

| | Echimidine | Integer-rimine | Intermedine | Senecionine | Seneci-phylline | Senkirkine |
|-------------------------------|------------|----------------|-------------|-------------|-----------------|------------|
| x_{PT} [µg/kg] | 70.8 | 54.4 | 43.9 | 29.1 | 12.9 | 62.1 |
| $u(x_{PT})$ [µg/kg] | 6.3 | 3.8 | 4.4 | 1.8 | 1.0 | 9.6 |
| \bar{x}_{rob} [µg/kg] | 68.4 | | 30.9 | 23.4 | 9.7 | 52.3 |
| σ_{PT} [µg/kg] | 15.6 | 12.0 | 9.7 | 6.4 | 2.8 | 13.7 |
| $u(x_{PT}) > 0.3 \sigma_{PT}$ | YES | YES | YES | NO | YES | YES |
| σ'_{PT} [µg/kg] | 16.8 | 12.6 | 10.6 | 6.6 | 3.0 | 16.7 |
| Number reported | 20 | 4 | 20 | 22 | 22 | 23 |
| "<" or ">" | 1 | 2 | 1 | 1 | 3 | 1 |
| $z' \leq 2 $ | 17 | 2 | 17 | 9 | 13 | 21 |
| $ 2 < z' \leq 3 $ | 1 | 0 | 2 | 3 | 3 | 0 |
| $z' > 3 $ | 1 | 0 | 0 | 9 | 3 | 1 |

Table 3 Reported values and calculated z'-scores for the PAs content in test item HT

| Lab Code | Integerrimine | | Retrorsine | | Senecionine | | Senecivernine | |
|----------|-----------------------------------|----------|---------------------|----------|---------------------|----------|---------------------|----------|
| | reported [µg/kg] | z'-score | reported [µg/kg] | z'-score | reported [µg/kg] | z'-score | reported [µg/kg] | z'-score |
| 1 | | | | | | | | |
| 2 | < 10 | | < 10 | | 75 | -0.3 | < 5 | |
| 3 | retronecine equivalents: 78 µg/kg | | | | | | | |
| 4 | | | 17.2 | -0.6 | 59.4 | -1.1 | 8.16 | 1.5 |
| 5 | | | 15.7 | -0.9 | 47.7 | -1.7 | 12 | 3.7 |
| 6 | < 2 | | 16.58 | -0.7 | 78.08 | -0.1 | 7.84 | 1.3 |
| 7 | | | 3.7 | -3.5 | 47.9 | -1.7 | | |
| 8 | | | 11 | -1.9 | 58.14 | -1.2 | 24.4 | 11.0 |
| 9 | | | 16.92 | -0.6 | 68.21 | -0.6 | 6.86 | 0.7 |
| 10 | | | < 10 | | 130 | 2.6 | 14.3 | 5.1 |
| 11 | | | 7.06 | -2.8 | 47.54 | -1.7 | 47.54 | 24.6 |
| 12 | | | 10.11 | -2.1 | 62.62 | -0.9 | | |
| 13 | 21.5 | -1.3 | 11.5 | -1.8 | 38.5 | -2.2 | 8.47 | 1.6 |
| 14 | | | < 2.5 | | 64 | -0.9 | | |
| 15 | | | 20.3 | 0.1 | 66.9 | -0.7 | < 5 | |
| 16 | | | < 2.5 | | 26.6 | -2.8 | | |
| 17 | | | < 3 | | 59.2 | -1.1 | | |
| 18 | | | | | 120 | 2.1 | 5.9 | 0.1 |
| 19 | | | | | | | | |
| 20 | | | 19.7 | 0.0 | 111 | 1.6 | 36 | 17.8 |
| 21 | | | | | 30 | -2.6 | | |
| 22 | | | < 10 | | 93 | 0.7 | 39 | 19.6 |
| 23 | | | < 20 | | < 30 | | | |
| 24 | | | 35.9 | 3.5 | 86.4 | 0.3 | 21.2 | 9.1 |
| 25 | | | | | 92.91 | 0.7 | | |
| 26 | 36.4 | 0.8 | 14.8 | -1.1 | 68.3 | -0.6 | 4.9 | -0.5 |

For this test item the measurements of the PAs Echimidine, Intermedine, and Senkirkine showed very satisfactory results with 85 % (17 of 20), 85 % (17 of 20), and 91 % (21 of 23), respectively, of the z'-scores below or equal to |2|. Analogous to item HT, only four laboratories reported for Integerrimine and, contrary to test item HT, the results for Senecionine in this item were very unsatisfactory.

Of the 22 z'-scores calculated for Senecionine nine (41 %) were larger than 3. This is all the more surprising since five of those nine laboratories (z'-scores 3.9 to 14.2) scored well for Senecionine in item HT (z'-scores -1.7 to 0.7). Additionally, seven laboratories (1, 5, 8, 10, 11, 20, 22) reported to have found Senecivernine in item HO and, of these, four scored well in this item for Senecionine with z'-scores between -1.0 and 0.9. The results are listed in Table 4 and visualised in Annex 6: Graphs for Test item HO. The results for Seneciphylline are also not convincing with only 13 out of 22 (59 %) z'-scores below or equal to |2|.

Laboratories 14 (83 µg/kg) and 19 (30 µg/kg) have reported to have found Lycopsamine and did not report results for Intermedine. Intermedine / Lycopsamine / Indicine form another group of

isomeric PAs (C₁₅H₂₅NO₅) which are difficult to separate. If laboratory 19 which indicated in the questionnaire that it reported Intermedine / Lycopsamine / Indicine as sum misidentified Intermedine as Lycopsamine then the reported value would have earned a z'-score of -1.3. Trichodesmine was reported to have been found by laboratory 12 (49 µg/kg). No other laboratory reported any other PA above its respective LOQ. Laboratory 3 reported the sum of all PAs as "retronecine equivalents" at 45 µg/kg.

Correct determination of the individual members of the isomeric groups "Integerrimine / Senecionine / Senecivernine" and "Intermedine / Lycopsamine / Indicine" seems to be an issue which needs attention.

Table 4 Reported values and calculated z'-scores for the PAs content in test item HO

| Lab Code | Echimidine | | Integerrimine | | Intermedine | | Senecionine | | Seneciphylline | | Senkirkine | |
|----------|-----------------------------------|----------|------------------|----------|------------------|----------|------------------|----------|------------------|----------|------------------|----------|
| | reported [µg/kg] | z'-score | reported [µg/kg] | z'-score | reported [µg/kg] | z'-score | reported [µg/kg] | z'-score | reported [µg/kg] | z'-score | reported [µg/kg] | z'-score |
| 1 | 57.6 | -0.8 | | | 24.7 | -1.8 | 74.7 | 6.9 | 9.1 | -1.3 | 74 | 0.7 |
| 2 | | | | | | | | | | | | |
| 3 | retronecine equivalents: 45 µg/kg | | | | | | | | | | | |
| 4 | 61.8 | -0.5 | | | 34.4 | -0.9 | 22.7 | -1.0 | 10 | -1.0 | 47.9 | -0.9 |
| 5 | 69.3 | -0.1 | | | 27.6 | -1.5 | 22.8 | -1.0 | 10.4 | -0.8 | 54.4 | -0.5 |
| 6 | 67.87 | -0.2 | < 0.5 | | 39.75 | -0.4 | 21.99 | -1.1 | 12.37 | -0.2 | 57.89 | -0.3 |
| 7 | 34 | -2.2 | | | 32.3 | -1.1 | 83 | 8.1 | 25.3 | 4.1 | 37.2 | -1.5 |
| 8 | 81.9 | 0.7 | | | 21.4 | -2.1 | 27.6 | -0.2 | 31.7 | 6.3 | 60.4 | -0.1 |
| 9 | 72.15 | 0.1 | | | 38.24 | -0.5 | 27.13 | -0.3 | 7.61 | -1.8 | 37.82 | -1.5 |
| 10 | 77 | 0.4 | | | 29.4 | -1.4 | 50 | 3.1 | 21.1 | 2.7 | 58.2 | -0.2 |
| 11 | 45.92 | -1.5 | | | 16.84 | -2.6 | 11.07 | -2.7 | 19.71 | 2.3 | 42.05 | -1.2 |
| 12 | 69.45 | -0.1 | | | 29.18 | -1.4 | 58.84 | 4.5 | 0.34 | -4.2 | 9.17 | -3.2 |
| 13 | 82.5 | 0.7 | 37.4 | -1.4 | 24.8 | -1.8 | 13.6 | -2.3 | 5.78 | -2.4 | 48.9 | -0.8 |
| 14 | | | | | | | 120 | 13.7 | 18 | 1.7 | 69 | 0.4 |
| 15 | 74.6 | 0.2 | | | 30.2 | -1.3 | 19.5 | -1.4 | 10.8 | -0.7 | 49.9 | -0.7 |
| 16 | 62.8 | -0.5 | | | 33 | -1.0 | 14.5 | -2.2 | 9.9 | -1.0 | 40.4 | -1.3 |
| 17 | | | | | 48 | 0.4 | 116.33 | 13.1 | < 3 | | 54.73 | -0.4 |
| 18 | 57 | -0.8 | | | | | 52 | 3.4 | 7.9 | -1.7 | 55 | -0.4 |
| 19 | 122 | 3.1 | < 1 | | | | 52 | 3.4 | < 1 | | 53 | -0.5 |
| 20 | 96.4 | 1.5 | | | 36.2 | -0.7 | 28.5 | -0.1 | 13.4 | 0.2 | 54.7 | -0.4 |
| 21 | 62 | -0.5 | | | 31 | -1.2 | | | | | 44 | -1.1 |
| 22 | 56 | -0.9 | | | 31 | -1.2 | 35 | 0.9 | 9 | -1.3 | 49 | -0.8 |
| 23 | > 10 | | | | > 10 | | > 10 | | > 10 | | > 10 | |
| 24 | | | | | | | | | | | | |
| 25 | | | | | 44.49 | 0.1 | 54.31 | 3.8 | 16.37 | 1.2 | 71.75 | 0.6 |
| 26 | 74 | 0.2 | 39.4 | -1.2 | 57 | 1.2 | 16.2 | -1.9 | 10.2 | -0.9 | 49.2 | -0.8 |

6.3. "Smaller than" / "Larger than" values and uncertainties

For the PAs listed in Table 1 and Table 2 16 results were reported as "smaller than". This declaration was only correct for three of the 16. For all others the actual mass fraction was

above the "smaller than" value. Five results were reported as "larger than" and all those declarations were correct.

Only for a minority (72 out of 173, i.e. 42 %) of results an associated measurement uncertainty was reported. Of those 72 expanded uncertainties 37 were under-estimated, meaning the uncertainty ranges around the reported value did not include the reference value. Another five reported uncertainty ranges included the value zero meaning the reported value within its uncertainty is not different from zero at the selected confidence range. In seven cases the estimate of the measurement uncertainty was overly conservative.

All affected laboratories are advised to check their LOD/LOQ and uncertainty estimations.

6.4. The questionnaire

All participating laboratories were asked to fill in a questionnaire addressing their accreditation and experience, preparation conditions for the two test items, chromatographic separation conditions, detection conditions, calibration approach, and any comments they felt like making. Each laboratory reporting results also filled in the questionnaire. For selected questions the correlation of answers to the z'-score of Senecionine results in test item HO or Senecivernine results in test item HT was analysed by single-factor analysis of variance (ANOVA) after variance stabilizing transformation of the response variable.

The vast majority of the laboratories (85 %) had more than 12 months of experience with PA analysis. Two laboratories reported to have less than 3 months and another two between 3 and 12 months of experience. Accreditation to ISO 17025 for PA analysis was held by 15 laboratories while 11 were not accredited. A question begging to be asked is whether holding accreditation has a positive effect on the reported results. Looking at the z'-scores of Senecionine results no difference in performance between accredited and non-accredited laboratories could be detected with ANOVA. The same is true for Senecivernine results. The data is depicted in Figure 6- 11 in Annex 6: Graphs.

The information provided about test item preparation showed that the laboratories used a median test portion size of 5 g for HO and 2 g for HT. This was then extracted with a median volume of 20 mL and 40 mL, for HO and HT, respectively. For both test items the median extraction time was 30 min. The composition of the extraction solvents for test item HO was either acidic aqueous (16), acidic aqueous/organic (2), neutral aqueous/organic (3), pure organic (1), or in one other case basic aqueous/organic. Three laboratories chose to not make a statement. For test item HT either acidic aqueous (19), acidic aqueous/organic (2), or neutral aqueous/organic (2) conditions were employed. Also here, three laboratories did not make statements. Neither for Senecionine results in HO nor for Senecivernine results in HT a significant effect of the extraction solvent composition on the z'-scores was detected by ANOVA.

In the next group of questions about chromatographic separation conditions the 26 laboratories all stated to have used liquid chromatography for separation of the PAs/PANOs. Methanol was used as organic mobile phase modifier by 19 and acetonitrile by five laboratories. Two laboratories did not provide information. Only in two cases the mobile phase did not contain any additives. Formate (formic acid with or w/o ammonium formate) was used by 17 participants, acetate (acetic acid with ammonium acetate) by 4, and one laboratory employed a basic mobile

phase with ammonium carbonate. Again, two laboratories did not provide information. Neither the organic modifier nor the additive in the mobile phase had a significant effect on the results of Senecionine in HO or Senecivernine in HT.

A hypothesis was raised above that the unsatisfactory result for Senecionine in HO was caused by insufficient chromatographic separation within the isomeric group "Integerrimine / Senecionine / Senecivernine". Ten laboratories reported to have used analytical columns with fully porous particles of larger than 2 μm size. The same number of laboratories used sub-2 μm fully porous particle columns. Since the sub-2 μm columns are capable of higher separation efficiencies it was tested whether that had a significant effect on the scores. But neither for Senecionine results in HO nor for the ones of Senecivernine in HT such an effect could be shown.

All laboratories have used mass spectrometry (MS) for detection with 20 laboratories using MS/MS in selected/multiple reaction monitoring (SRM/MRM), one using high resolution MS/MS, and three using selected ion monitoring (SIM). Two made no statement. The question how the PAs/PANOs were determined was answered by 24 laboratories with individual PAs and PANOs as such. One laboratory reduced all PAs/PANOs to a common analyte (retronecine) and derivatised with phthalic acid before measurement. And only one laboratory (Lab 19) stated that the sums of isomers were reported.

Questions about the calibration approach were not included in the initial questionnaire. Therefore, all 26 reporting laboratories received an email with these questions on 09.11.2017. All except one laboratory replied. For the test item HO 24 laboratories reported results and, of those, 12 used matrix-matched external standard calibration, six standard addition, three neat solvent external standard calibration, one isotope dilution, one "neat solvent external standard/standard addition", and one "prespiked matrix matched calibration". The same number (24) of laboratories reported at least one result for test item HT. Of those, nine used standard addition, six matrix-matched external standard calibration, five neat solvent external standard calibration, one isotope dilution, one "neat solvent external standard/standard addition", and one "prespiked matrix matched calibration". One of these laboratories did not reply to the inquiry. The influence of the calibration approach on the z' -scores for Senecionine in HO and Senecivernine in HT was tested with ANOVA and no difference could be found.

Altogether, 21 laboratories reported results for both test items and only five of them reported all their results with satisfactory z' -scores. Laboratory 26 using standard addition reported ten values with z' -scores between -1.1 and 2.0. Laboratory 6 reported eight values with scores between -1.1 and 1.3 and two "smaller than" results. Laboratory 4 reported eight values with scores between -1.1 and 1.5. Both laboratories used matrix-matched external standard calibration for test item HO and standard addition for HT. Laboratory 9 reported eight values between -1.8 and 0.7 using neat solvent external standard calibration. Finally, laboratory 15 reported seven values with scores between -1.5 and 0.2 using standard addition and one "smaller than" result.

In the comments section of the questionnaire laboratory 11 declared that it reported the sum of Senecionine and Senecivernine. It reported the same value for both in HT. That led to a satisfactory z' -score for Senecionine and an unsatisfactory one for Senecivernine. That Indicine and Intermedine could not be separated on a C18 column was stated by laboratory 8. Laboratory

14 commented that co-occurrence of the isomeric PAs might occur and affect their results. Co-elution of Lycopsamine and Indicine was remarked by laboratory 16.

7. Conclusions

This proficiency testing scheme was executed to assess the capabilities of laboratories to determine PAs/PANOs in herbal infusions and honey. 29 laboratories from nine EU Member States plus Singapore registered for this PT scheme and were sent two test items for analysis. By the dead line of 24.10.2017 26 laboratories had reported back results and filled in a questionnaire.

Only three laboratories out of the 26 were able to test for all 10 measurands and only one reported all 10 values with z'-scores smaller or equal to $|2|$. Overall only five laboratories obtained satisfactory z'-scores ($\leq |2|$) for all their reported results. The outcome of this scheme points to a possible problem by many laboratories to properly determine the individual members of the two isomeric groups "Integerrimine / Senecionine / Senecivernine" and "Intermedine / Lycopsamine / Indicine".

For test item HO, spiked with six PAs/PANOs, the measurements of the PAs Echimidine, Intermedine, and Senkirkine showed very satisfactory results with 85 % (17 of 20), 85 % (17 of 20), and 91 % (21 of 23), respectively, of the z'-scores below or equal to $|2|$. For the four PAs present in test item HT a satisfying outcome can only be attested for the analysis of Senecionine. It was present at the highest level and 17 out of 23 (74 %) reported results showed z'-scores $\leq |2|$.

8. Acknowledgements

We want to thank the participating laboratories (see Annex 8: Participating laboratories, p.44) for their efforts in reporting the results on time and filling in the questionnaire. Our thanks also go to C. Mischke of the EURL for Mycotoxins for his support in test material production, and Vasiliki Exarchou for the qNMR measurements.

References:

- 1 EFSA (European Food Safety Authority), 2016. Dietary exposure assessment to pyrrolizidine alkaloids in the European population. *EFSA Journal* 2016;14:4572, 50 pp. doi:10.2903/j.efsa.2016.4572
- 2 ISO/IEC 17043 "Conformity assessment – General requirements for proficiency testing, International Organisation for Standardization, Geneva, Switzerland, 2010.
- 3 ISO 13528:2015 - Statistical methods for use in proficiency testing by interlaboratory comparison, International Organisation for Standardization, Geneva, Switzerland, 2015.
- 4 Hauswaldt, A.-L.; Rienitz, O.; Jährling, R.; Fischer, N.; Schiel, D.; Labarraque, G.; Magnusson, B. Uncertainty of standard addition experiments: a novel approach to include the uncertainty associated with the standard in the model equation. *Accreditation and Quality Assurance* 2012, 17, 129-138, 10.1007/s00769-011-0827-5.
- 5 Thompson, M.; Ellison, S.L.R.; Wood, R. The International Harmonized Protocol for the Proficiency Testing of Analytical Chemistry Laboratories (IUPAC Technical Report). *Pure Appl. Chem.* **2006**, 78, 145-196.

Annex 1: Instructions for PT scheme



EUROPEAN COMMISSION
 DIRECTORATE-GENERAL
 JOINT RESEARCH CENTRE
 Directorate F - Health, Consumers & Reference Materials (Geel/Ispra)
 Food & Feed Compliance

Instructions for PT scheme

"Pyrrolizidine alkaloids contamination in herbal tea and honey"

«AddressBlock»

«GreetingLine»

please read the following instructions carefully.

This PT scheme is to assess your laboratory's ability to determine pyrrolizidine alkaloids (PAs) in herbal tea and/or honey. The following PAs may or may not be present in the two test materials: Echimidine, Erucifoline, Europine, Indicine, Integerrimine, Intermedine, Jacobine, Lasiocarpine, Lycopsamine, Monocrotaline, Retrorsine, Senecionine, Seneciphylline, Senecivernine, Senkirkine, Trichodesmine.

It is your choice whether you measure just one or both of the two test materials and which of the 16 PAs you determine. You may report a sum parameter combining all PAs or results for the individual PAs. If results for individual PAs are reported you must report the sum of the amine and the N-oxide of the respective PA. This sum can either be calculated from individual values for the amine and N-oxide or it may be the measured value of the total amine after reduction of the N-oxide.

Please check the integrity of the test material containers and use the included material receipt form to report back to us your findings as soon as possible. Until commencement of the analysis the test materials should be stored at deep freeze temperatures (-15 - -25 °C). Before opening the containers they should be given sufficient time to equilibrate to room temperature. It is the responsibility of the laboratory to withdraw a representative test portion from the test material containers. Therefore, to minimize sub sampling uncertainties, the content of the container should be mixed well before withdrawal and the test portion be withdrawn in multiple smaller increments.

Scoring of the reported results will be done acc. to ISO 13528:2015 and be based on an independently determined reference value (x_{PT}) for each PA and a target standard deviation σ_{PT} of $0.22 * x_{PT}$ acc. to the following equation:

$$z = \frac{(x_i - x_{PT})}{\sigma_{PT}}$$

with x_i the reported result of an individual PA from a laboratory. No scoring will be provided for reported sums of total PAs.

The reporting deadline is **20. Oct 2017**. After this date the reporting interface will be closed.

Open the following link (<https://web.jrc.ec.europa.eu/ilcReportingWeb/>) in a web browser (preferably use latest version Firefox or Internet Explorer) to enter your results.

Your personal password key is: «Part_key»

Enter this into the "Password key" field.

When entering your results use a "." (dot) as decimal divider. Do not use a "," (comma). To enter a result below LOD into the results table select the "<" (smaller as) sign from the drop down menu and enter the numerical value of the LOD into the results field. Then choose "OTHER, Please specify" in the Technique column and enter "LOD".

For example, a measurement of Erucifoline with a LOD of 3 µg/kg would be entered as is depicted below:

| Measurand | Measurement | Reference Date | Result | Unit | Uncert. value | Coverage Faktor k | Technique | Clear |
|------------------|----------------------|----------------|--------|-------|---------------|-------------------|------------------------------|-------|
| ! Sum of all PAs | Concentration[µg/kg] | Value | = | µg/kg | | | No technique | |
| Echimidine | Concentration[µg/kg] | Value | = | µg/kg | | | No technique | |
| Erucifoline | Concentration[µg/kg] | Value | < | µg/kg | | | OTHER, Please specify LOD | |
| Europine | Concentration[µg/kg] | Value | = | µg/kg | | | No technique | |

For entering a value below LOQ proceed accordingly. Leave the fields for PAs for which no measurements were performed empty. If you determined the sum of all PAs enter your result into the top row and leave all others empty. Results entry does not need to happen in one session. Entry can always be interrupted. Save what you have entered and come back later to continue. Only once you click the "Submit my results" button on the reporting home page the entries will be submitted.

If possible also provide an estimate of the measurement uncertainty and the coverage factor used for any reported results larger than LOQ. The annex provides further information on how to use the reporting interface.

There is also a questionnaire about your analysis that we ask you to fill in at the time of reporting.

To prevent falsification of results you must keep the password key confidential. Any attempt of collusion is contrary to professional scientific conduct and serves only to nullify the benefits of proficiency tests to costumers, accreditation bodies and analysts alike.

Good luck with your analysis.

Andreas Breidbach
PT coordinator

Annex 2: Proficiency testing materials receipt form



EUROPEAN COMMISSION

DIRECTORATE-GENERAL
JOINT RESEARCH CENTREDirectorate F - Health, Consumers and Reference Materials
Food & Feed Compliance

Geel, 12 December 2017

PROFICIENCY TESTING MATERIALS RECEIPT FORM**Name:** «Title» «Firstname» «Surname»**Institute:** «Organisation»**NOTE: STORE MATERIAL IN A FREEZER AT -15 - -25 °C!**

Please ensure that the items listed below have been received undamaged, and then check the relevant statement:

| | |
|--|----------|
| Date of receipt | |
| All items have been received undamaged | YES / NO |
| <i>If NO, please list damaged items:</i> | |

Contents of the parcel:

- a) 1 container with artificially contaminated honey (HO)
- b) 1 container with naturally contaminated herbal tea (HT)
- c) Instructions with information about storage and reporting
- d) this material receipt form

Please scan in the completed form and e-mail to:

JRC-EURL-MYCOTOX@ec.europa.eu

Your Signature / Stamp here:

Annex 3: Homogeneity testing

For homogeneity testing of the HT test material 10 test units were randomly selected out of the 96 test units that were created. Of each of the 10 test units two 1 g test portions were weighed out for extraction and clean-up. Care was taken to ascertain that each test portion was representative of its respective test unit.

For homogeneity testing of the HO test material eight test units were selected from the 73 test units that were created. Since the test units were created from a more or less continuous material stream a random unit was selected from the first 10 units. Then every 10th following unit was selected until the end. In this particular case the units 3, 13, 23, 33, 43, 53, 63, and 73 were selected. Of each of the selected units 5 g were weighed into a new tube. Then 2.5 g of water were added to reduce viscosity of the honey. After thorough mixing two 3 g test portions, representing 2 g test material, were weighed into new tubes for clean-up.

The test portions were extracted/cleaned-up with SPE acc. to an internal protocol and measured in randomized order with a LC-HRMS. The obtained signals, without conversion to concentration units, were evaluated for sufficient homogeneity acc. to ISO 13528 Annex B [3]. Table 1- 1 shows the results for test material HT, Table 1- 2 the results for material HO.

Table 1- 1: Results of the homogeneity test for test material HT

| Analyte | c | s_s^2 | s_w^2 |
|---------------|---------|---------|---------|
| Integerrimine | 0.178 | 0.0133 | 0.0905 |
| Retrorsine | 0.0142 | 0 | 0.00742 |
| Senecionine | 0.93 | 0.059 | 0.55 |
| Senecivernine | 0.00927 | 0.00234 | 0.00391 |

c - critical value, s_s^2 - between-sample variance, s_w^2 - within-sample variance

Table 1- 2: Results of the homogeneity test for test material HO

| Analyte | c | s_s^2 | s_w^2 |
|----------------|--------|----------|----------|
| Echimidine | 2.74 | 0 | 0.310 |
| Integerrimine | 1.08 | 0 | 0.0255 |
| Intermedine | 0.558 | 0.000472 | 0.0207 |
| Senecionine | 0.272 | 0.00101 | 0.00148 |
| Seneciphylline | 0.0493 | 0.000464 | 0.000899 |
| Senkirkine | 2.20 | 0.0221 | 0.0909 |

c - critical value, s_s^2 - between-sample variance, s_w^2 - within-sample variance

For all analytes in both test materials the between-sample variance was smaller than the critical value indicating sufficient homogeneity.

Annex 4: Stability testing

The stability of the analytes in the two test materials was determined through an isochronous study. For this purpose randomly selected test units were stored frozen (-15 to -25 °C). Eight weeks before the measurement a set of test units was transferred to and stored at 40 and 4 °C, respectively. Four weeks before the measurement another set of test units was transferred to and stored at the same temperatures. At the day of measurement a set of test units stored frozen, the test units stored for 4 and 8 weeks at 40 °C, and the test units stored for 4 and 8 weeks at 4 °C were prepared and measured as described for the homogeneity testing. The test units stored frozen represented time point zero.

Based on the statistical model with y the measured mass fraction at time point x (in weeks)

$$y = \beta_0 + \beta_1 x$$

the coefficients β_0 and β_1 were computed via least square regression. The slope β_1 represents the analyte change per time unit. If the 95 % confidence range around β_1 includes zero there is no statistically significant instability. Table 2- 1 shows that this is the case for all analytes in test material HT but Retrorsine stored at 40 °C. The temperature level "40 °C" was included to simulate the worst case scenario for transportation. All shipments arrived in the respective laboratory within 2 days. After reception laboratories were required to store the test materials at deep-freeze temperatures. For Retrorsine in HT a 2 day (0.29 weeks) exposure to transportation temperature would result in a analyte change of $-0.23 \text{ ng/g/week} \times 0.29 \text{ weeks} = -0.066 \text{ ng/g}$. This is smaller than $0.3 \sigma_p$ (1.2 ng/g) and has, therefore, no impact on the calculated z-scores. This is indicated in the right-most column of Table 2- 1.

Table 2- 1: Results of the stability study for test material HT

| Analyte | Temperature [°C] | β_1 [ng/g week] | 95 % CR [ng/g week] | $\beta_1 \times \max(t_T) > 0.3 \sigma_p$ |
|---------------|------------------|-----------------------|---------------------|---|
| Integerrimine | 4 | 0.029 | -0.75 to 0.81 | |
| | 40 | -0.22 | -0.68 to 0.23 | |
| Retrorsine | 4 | -0.16 | -0.62 to 0.31 | |
| | 40 | -0.23 | -0.46 to -0.0032 | NO |
| Senecionine | 4 | 0.41 | -1.2 to 2.0 | |
| | 40 | -0.77 | -1.9 to 0.35 | |
| Senecivernine | 4 | -0.092 | -0.29 to 0.11 | |
| | 40 | -0.068 | -0.23 to 0.10 | |

β_1 : contamination change per week, 95 % CR: 95 % confidence range of β_1 , $\beta_1 \times \max(t_T) > 0.3 \sigma_p$: $\beta_1 \times$ maximum possible exposure time to respective temperature is larger than $0.3 \times$ target standard deviation

Table 2- 2: Results of the stability study for test material HO

| Analyte | Temperature [°C] | β_1 [ng/g week] | 95 % CR [ng/g week] | $\beta_1 \times \max(t_T) > 0.3 \sigma_p$ |
|----------------|------------------|-----------------------|---------------------|---|
| Echimidine | 4 | -0.67 | -1.2 to -0.16 | NO |
| | 40 | -0.47 | -1.0 to 0.11 | |
| Integerrimine | 4 | -0.27 | -0.64 to 0.086 | |
| | 40 | 0.052 | -0.34 to 0.45 | |
| Intermedine | 4 | -0.36 | -0.69 to -0.028 | NO |
| | 40 | 0.015 | -0.29 to 0.32 | |
| Senecionine | 4 | -0.11 | -0.29 to 0.064 | |
| | 40 | -0.10 | -0.35 to 0.15 | |
| Seneciphylline | 4 | -0.010 | -0.14 to -0.060 | NO |
| | 40 | -1.2 | -1.6 to -0.80 | NO |
| Senkirkine | 4 | -0.026 | -0.33 to 0.27 | |
| | 40 | 0.073 | -0.38 to 0.53 | |

β_1 : contamination change per week, 95 % CR: 95 % confidence range of β_1 , $\beta_1 \times \max(t_T) > 0.3 \sigma_p$: $\beta_1 \times$ maximum possible exposure time to respective temperature is larger than $0.3 \times$ target standard deviation

Table 2- 2 lists the result for the test material HO. Here for Echimidine and Intermedine a statistically significant instability was detected at 4 °C storage but not at 40 °C. But even incorrect storage at 4 °C over the 7 week study period would not lead to a contamination change which would affect the z-score. For Seneciphylline instabilities at both tested temperatures were detected. Assuming a 2-day exposure to 40 °C during transportation and 7 weeks of storage at 4 °C an analyte change of 0.41 ng/g would result. The critical value $0.3 \times \sigma_p$ is 0.89 ng/g. All of this is indicated in the right-most column.

Annex 5: Determination of the assigned values and their uncertainties

This was done with gravimetric standard addition using an internal standard (ISTD) acc. to Hauswaldt et al. [4]. In brief, an amount of a mixed PAs solution was gravimetrically added to an aliquot of the test material. The mixed PAs solution was prepared from individual PA reference materials of which the purity was verified using HPLC, TGA-GC-MS, and qNMR. For the HT test material approx. 15 or 77 mg were added in duplicate to four aliquots. For the HO test material approx. 27, 56, or 80 mg were added to six aliquots in duplicate. Since the HO test material contained Seneciophylline-NO approx. 37 or 74 mg of a Seneciophylline-NO reference solution were added in duplicate to another four aliquots. For each test material two aliquots were prepared without addition of PAs. To all HT aliquots a constant amount of approx. 157 mg of an ISTD solution (Heliotrine) was added. To all HO aliquots approx. 192 mg ISTD were added. All additions were performed with a calibrated analytical balance with a readability of $d = 0.01$ mg.

Extraction and clean-up were performed the day after the additions to give the added reference materials sufficient time to equilibrate. To this end, H₂O/formic acid (95/5, v/v) and zinc were added to the aliquots. This created a reductive environment to reduce the N-oxides to the respective amines. After 60 min the zinc and other particulate matter were pelleted through centrifugation. The supernatant was adjusted to approx. pH 10 with conc. ammonia and an aliquot was cleaned-up with SPE. The dried eluate was reconstituted with H₂O/formic acid (999/1, v/v) and injected into a LC-HRMS. Almost baseline separation was achieved for the critical groups: intermedine/lycopsamine and senecivernine/integerrimine/Senecionine (see Figure 5- 1). The reductive extraction was chosen to reduce the complexity of the separation and to limit the number of reference materials necessary for this task.

Iteratively weighted least square regression was used to determine the intercept β_0 and the slope β_1 of the linear model:

$$R' \frac{m_y}{m_x} = \beta_0 + \beta_1 \frac{m_z}{m_x}$$

with

R' – ratio of peak area of analyte divided by peak area of ISTD;

m_y – mass of the ISTD solution;

m_x – mass of the test material aliquot;

m_z – mass of the mixed PAs solution.

With β_0 and β_1 known the mass fraction of the analyte w_x can be calculated:

$$w_x = \frac{\beta_0}{\beta_1} w_z$$

with w_z the mass fraction of the analyte in the mixed PAs solution.

Following the law of error propagation the uncertainty of w_x is

$$u(w_x) = \sqrt{\left(\frac{w_x}{w_z}\right)^2 u^2(w_z) - 2 \frac{w_z w_x}{\beta_1^2} u(\beta_0, \beta_1) - \left(\frac{w_x}{\beta_1}\right)^2 u^2(\beta_1) + \left(\frac{w_z}{\beta_1}\right)^2 u^2(\beta_0)}$$

with

$u^2(w_z)$ – squared standard uncertainty of the analyte mass fraction in the mixed PAs solution;

$u(\beta_0, \beta_1)$ – co-variance of the intercept β_0 and the slope β_1 from the regression;

$u^2(\beta_1)$ – squared standard error of the slope β_1 from the regression;

$u^2(\beta_0)$ – squared standard error of the intercept β_0 from the regression.

Acc. to ISO 13528:2015 [3] the assigned value for an analyte in a test material is calculated as

$$x_{PT} = w_x + \delta_{homo} + \delta_{stab} + \delta_{trans}$$

with

x_{PT} – the assigned value;

w_x – the reference value determined through standard addition;

δ_{homo} – error due to possible inhomogeneity (assumed as zero);

δ_{stab} – error due to possible instability during the proficiency test period (assumed as zero);

δ_{trans} – error due to possible instability under transport conditions.

Since those errors were negligible they were included into the respective standard uncertainties. The associated combined uncertainty is calculated as

$$u(x_{PT}) = \sqrt{u^2(w_x) + u^2(homo) + u^2(stab) + u^2(trans)}$$

with

$u(x_{PT})$ – combined uncertainty of the assigned value;

$u^2(w_x)$ – squared uncertainty of the mass fraction of the analyte;

$u^2(homo)$ – squared uncertainty of the homogeneity determination incl. possible error;

$u^2(stab)$ – squared uncertainty of the stability determination incl. possible error;

$u^2(trans)$ – squared uncertainty due to transport conditions incl. possible error.

For all analytes $u(stab)$ was calculated as

$$u(stab) = \frac{RSD_y}{\sum(x_i - \bar{x})} \times 7 \times w_x$$

with

RSD_y – relative standard deviation of all mass fractions determined during the stability study for a given analyte;

x_i – individual time points;

\bar{x} – average of all time points;

w_x – the reference value determined through standard addition;

the factor 7 represents the duration of the PT scheme in weeks.

For analytes for which a significant instability was detected the values for δ_{trans} and $u^2(trans)$ were calculated from the slope β_1 and its uncertainty $u(\beta_1)$, respectively, by multiplying with 2/7 (transportation time of 2 days expressed in weeks). The main contributors to the uncertainty budgets of the assigned values were the uncertainties due to the determination of β_0 and β_1 . Table 3- 1 lists the respective values for test material HT, and Table 3- 2 does the same for test material HO.

Table 3- 1: The assigned values and their uncertainties for the test material HT

| Analyte | x_{PT} [ng/g] | $u(w_x)$ [ng/g] | $u(homo)$ [ng/g] | $u(stab)$ [ng/g] | $u(trans)$ [ng/g] | $u(x_{PT})$ [ng/g] |
|----------------------|--------------------|--------------------|---------------------|---------------------|----------------------|-----------------------|
| Integerrimine | 30.6 | 2.8 | 0.11 | 0.14 | 0 | 2.8 |
| Retrorsine | 19.9 | 1.5 | 0 | 0.099 | 0.11 | 1.5 |
| Senecionine | 80.4 | 7.8 | 0.24 | 0.36 | 0 | 7.8 |
| Senecivernine | 5.7 | 1.2 | 0.048 | 0.037 | 0 | 1.2 |

Table 3- 2: The assigned values and their uncertainties for the test material HO

| Analyte | x_{PT} [ng/g] | $u(w_x)$ [ng/g] | $u(homo)$ [ng/g] | $u(stab)$ [ng/g] | $u(trans)$ [ng/g] | $u(x_{PT})$ [ng/g] |
|-----------------------|--------------------|--------------------|---------------------|---------------------|----------------------|-----------------------|
| Echimidine | 71.1 | 6.3 | 0 | 0.099 | 0.31 | 6.3 |
| Integerrimine | 54.4 | 3.8 | 0 | 0.076 | 0 | 3.8 |
| Intermedine | 44.1 | 4.4 | 0.022 | 0.088 | 0.16 | 4.4 |
| Senecionine | 29.1 | 1.8 | 0.032 | 0.036 | 0 | 1.8 |
| Seneciphylline | 13.4 | 1.0 | 0.021 | 0.29 | 0.51 | 1.2 |
| Senkirkine | 62.1 | 9.6 | 0.15 | 0.086 | 0 | 9.6 |

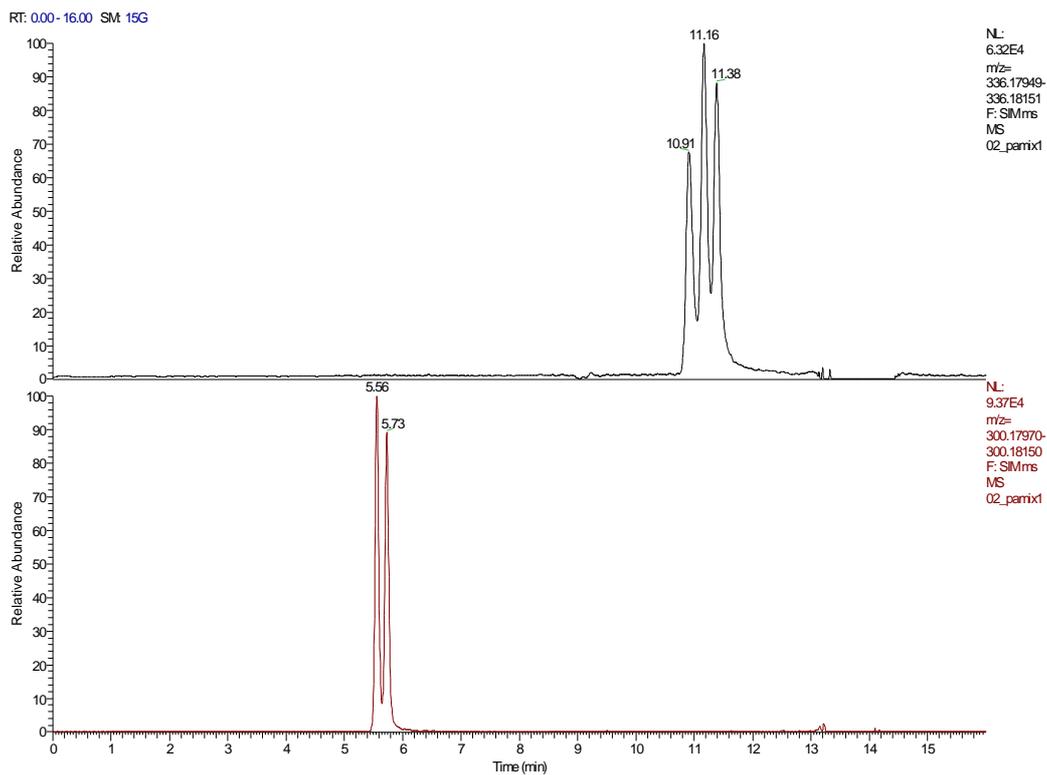


Figure 5- 1 Extracted ion chromatograms of a PA mixed reference solution; the top panel shows the separation of Senecivernine / Integerrimine / Senecionine (m/z 336.1805 \pm 3 ppm); the bottom panel shows the separation of Intermedine / Lycopsamine (m/z 300.1805 \pm 3 ppm); Indicine coelutes with Lycopsamine under these conditions; separation was afforded by a Phenomenex Kinetex C₁₈ 100x2.1mm, 1.7 μ m particle size.

Annex 6: Graphs

Test item HT

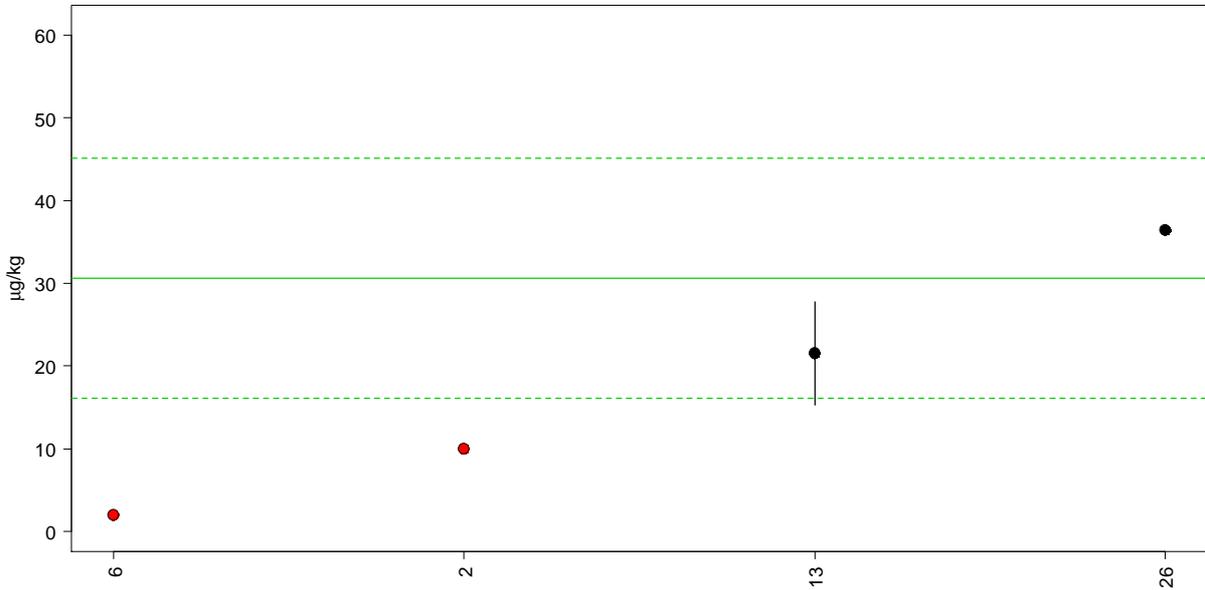


Figure 6- 1 Plot of results for Integerrimine from smallest to largest; black circles – reported results, red circles – reported "smaller than", black vertical line – reported measurement uncertainty, green solid line – assigned value, green broken line – two times adjusted target standard deviation σ'_{PT}

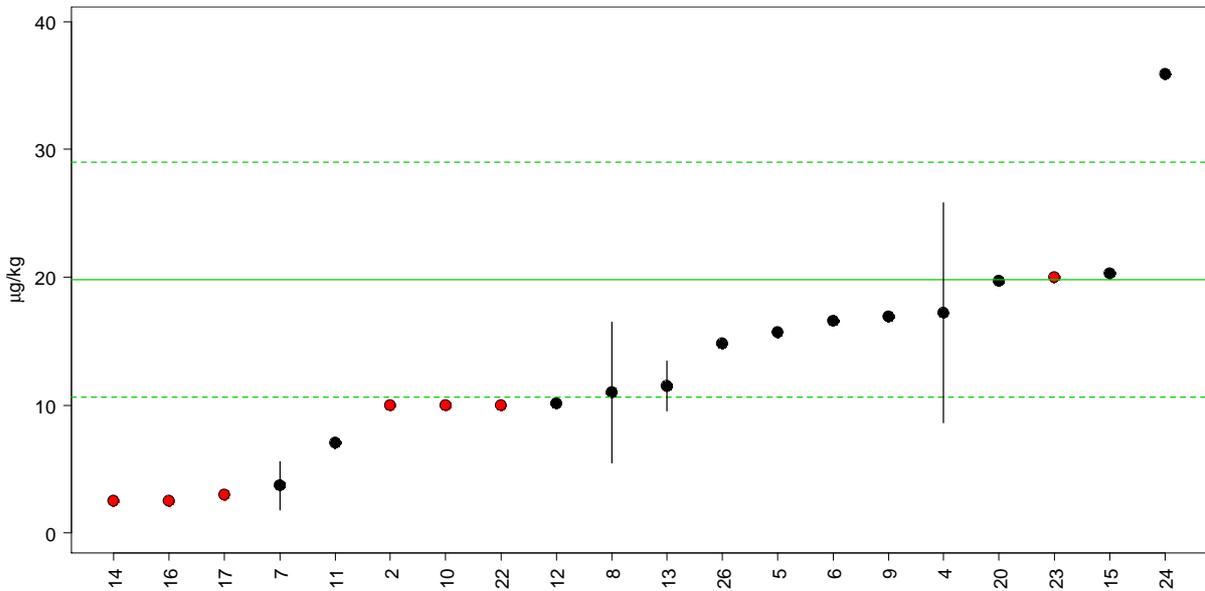


Figure 6- 2 Plot of results for Retrorsine from smallest to largest; black circles – reported results, red circles – reported "smaller than", black vertical line – reported measurement uncertainty, green solid line – assigned value, green broken line – two times adjusted target standard deviation σ'_{PT}

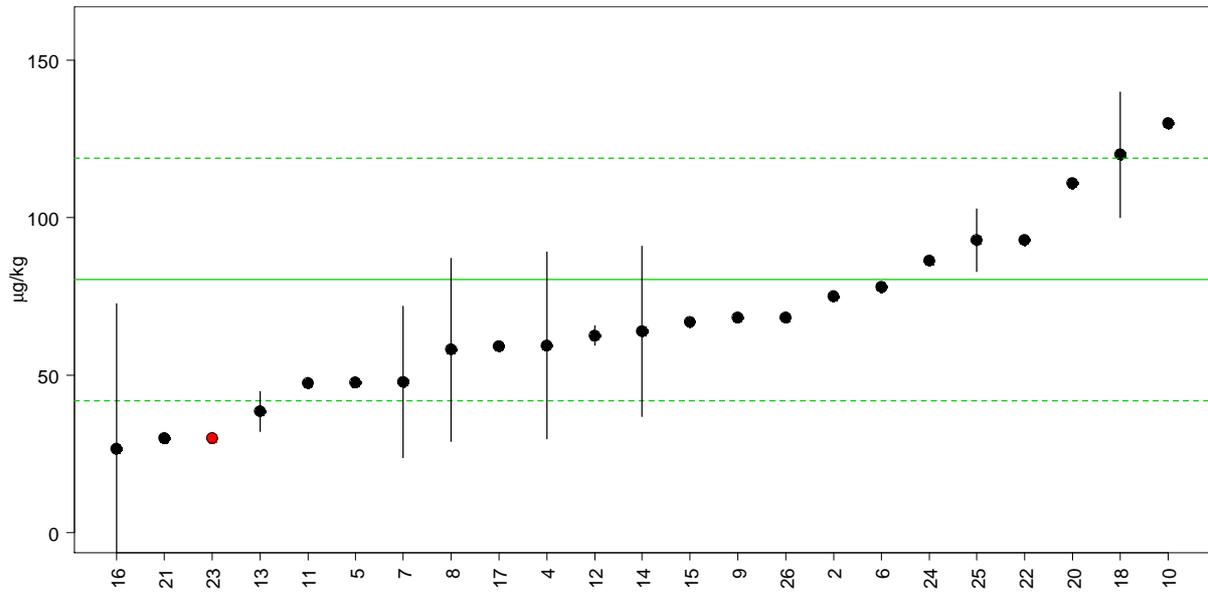


Figure 6- 3 Plot of results for Senecionine from smallest to largest; black circles – reported results, red circles – reported "smaller than", black vertical line – reported measurement uncertainty, green solid line – assigned value, green broken line – two times adjusted target standard deviation σ'_{PT}

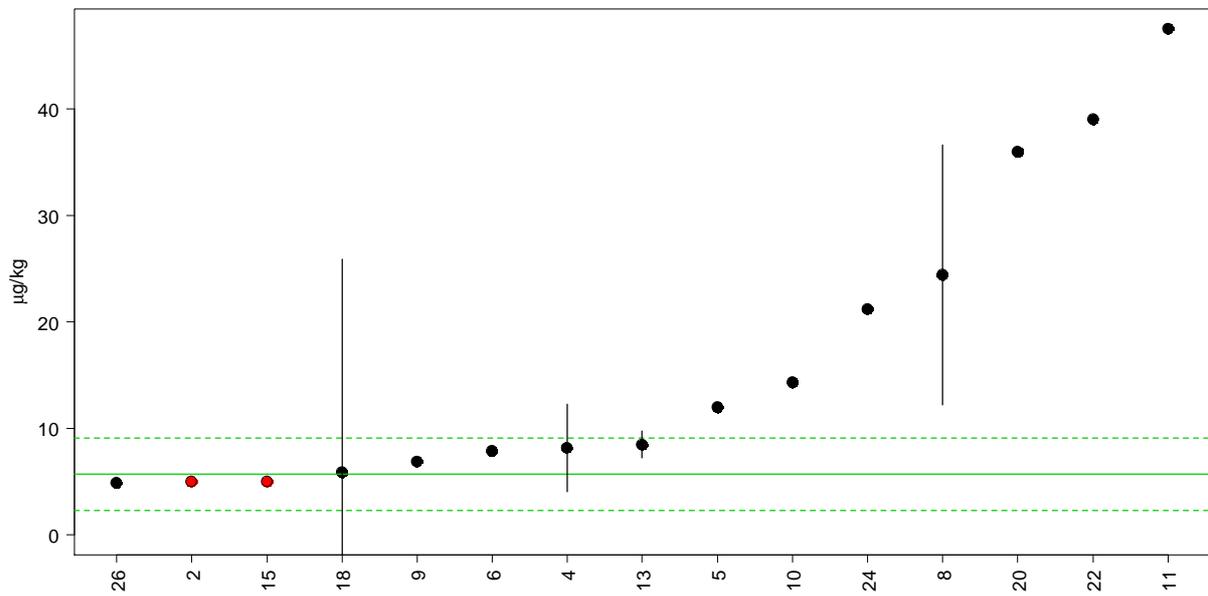


Figure 6- 4 Plot of results for Senecivernine from smallest to largest; black circles – reported results, red circles – reported "smaller than", black vertical line – reported measurement uncertainty, green solid line – assigned value, green broken line – two times adjusted target standard deviation σ'_{PT}

Test item H0

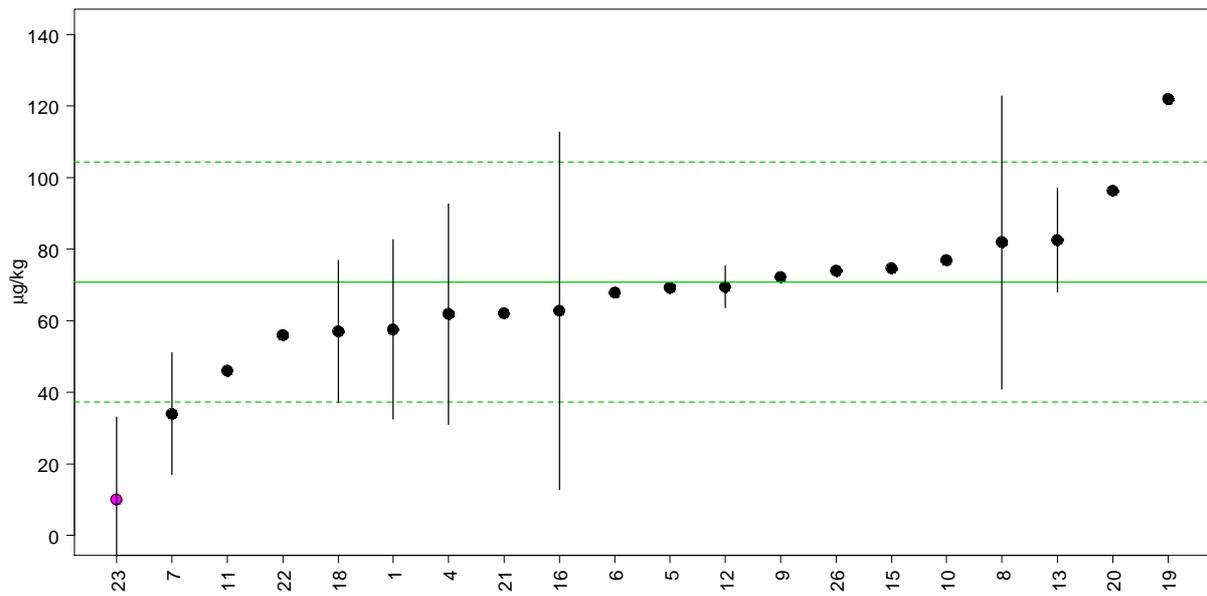


Figure 6- 5 Plot of results for Echimidine from smallest to largest; black circles – reported results, magenta colored circles – reported "larger than", black vertical line – reported measurement uncertainty, green solid line – assigned value, green broken line – two times adjusted target standard deviation σ'_{PT}

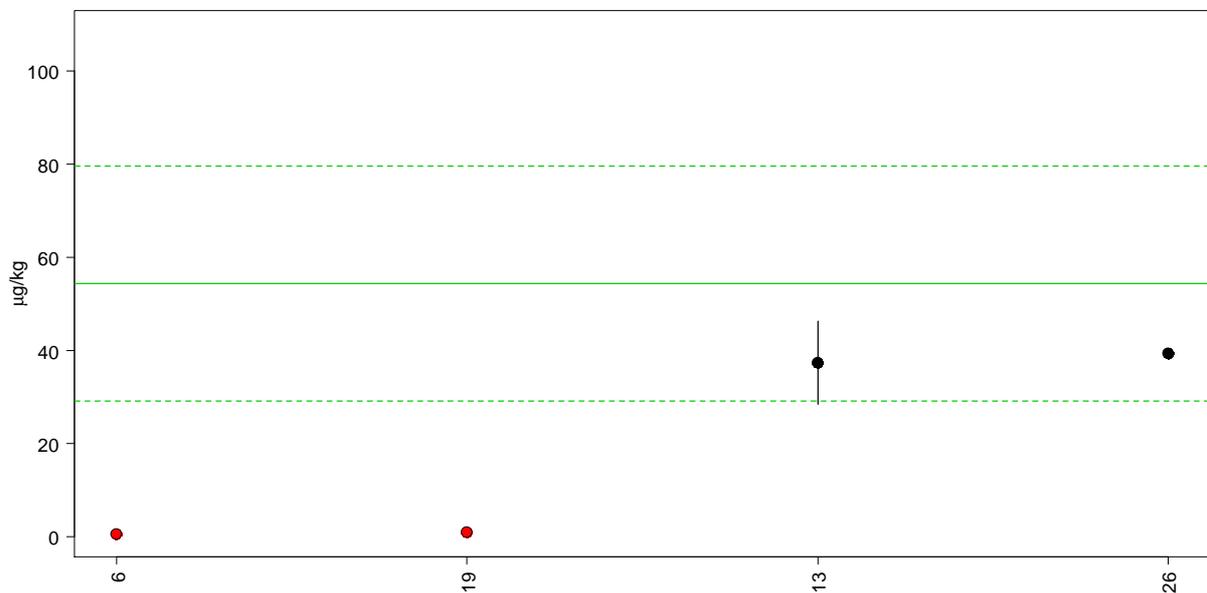


Figure 6- 6 Plot of results for Integerrimine from smallest to largest; black circles – reported results, red circles – reported "smaller than", black vertical line – reported measurement uncertainty, green solid line – assigned value, green broken line – two times adjusted target standard deviation σ'_{PT}

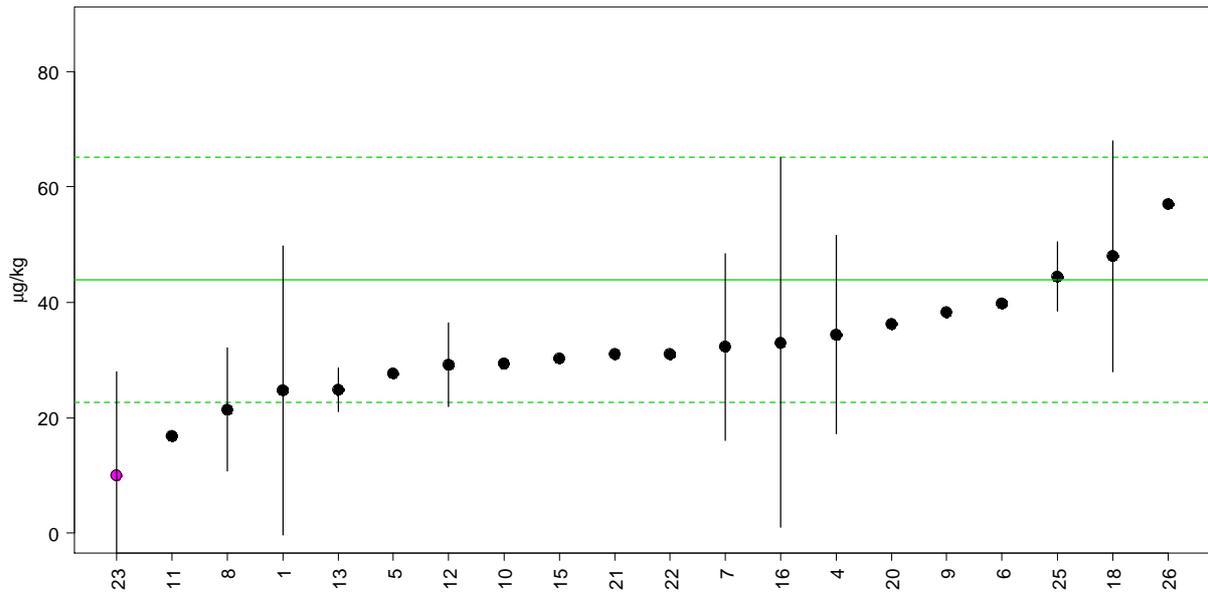


Figure 6- 7 Plot of results for Intermedine from smallest to largest; black circles – reported results, magenta colored circles – reported "larger than", black vertical line – reported measurement uncertainty, green solid line – assigned value, green broken line – two times adjusted target standard deviation σ'_{PT}

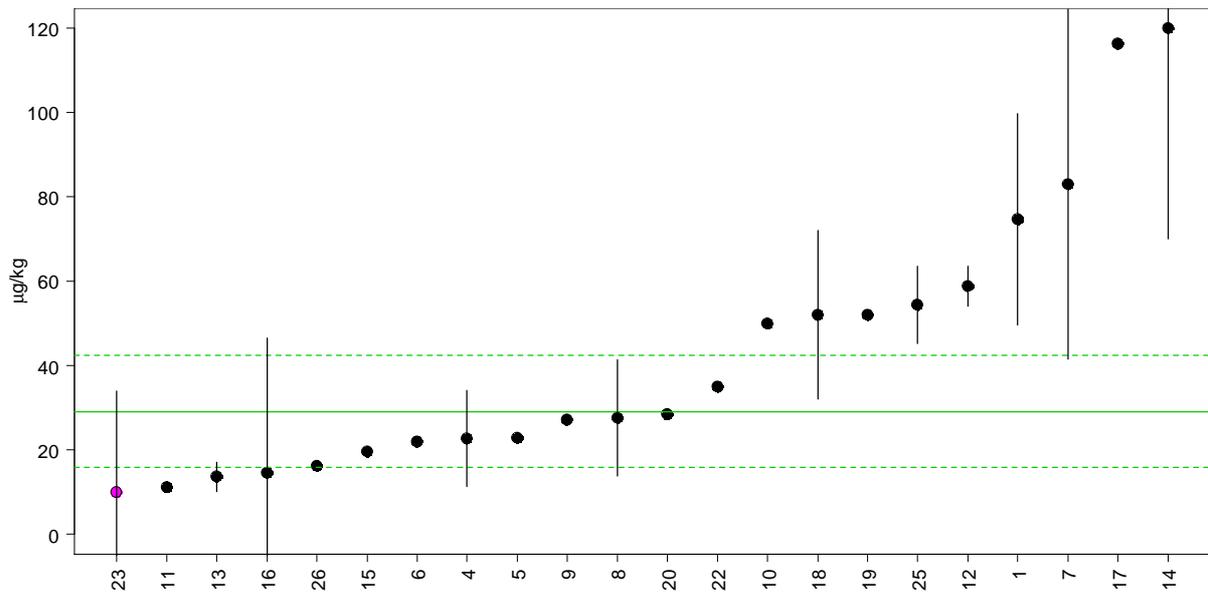


Figure 6- 8 Plot of results for Senecionine from smallest to largest; black circles – reported results, magenta colored circles – reported "larger than", black vertical line – reported measurement uncertainty, green solid line – assigned value, green broken line – two times adjusted target standard deviation σ'_{PT}

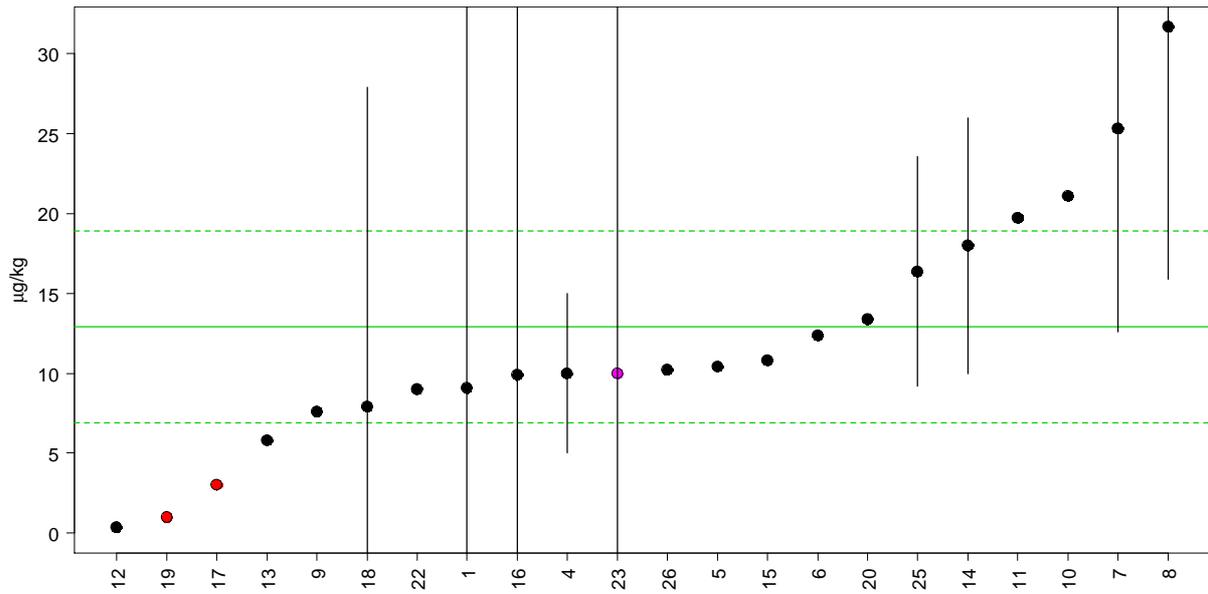


Figure 6- 9 Plot of results for Seneciphylline from smallest to largest; black circles – reported results, magenta colored circles – reported "larger than", red circles – reported "smaller than", black vertical line – reported measurement uncertainty, green solid line – assigned value, green broken line – two times adjusted target standard deviation σ'_{PT}

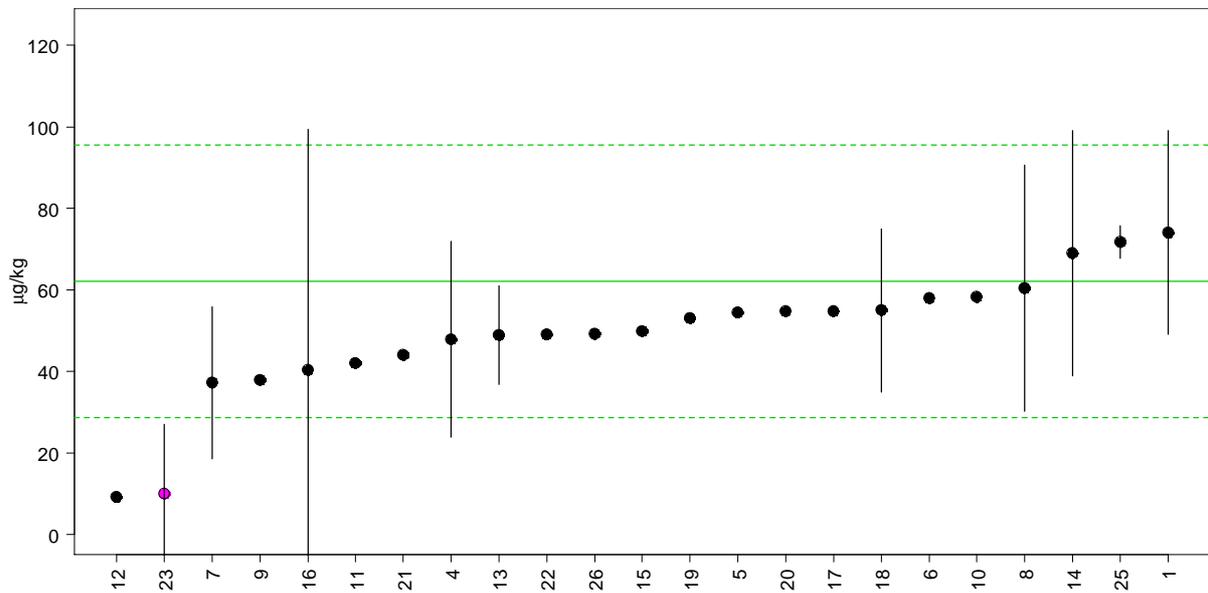


Figure 6- 10 Plot of results for Senkirikine from smallest to largest; black circles – reported results, magenta colored circles – reported "larger than", black vertical line – reported measurement uncertainty, green solid line – assigned value, green broken line – two times adjusted target standard deviation σ'_{PT}

Questionnaire data

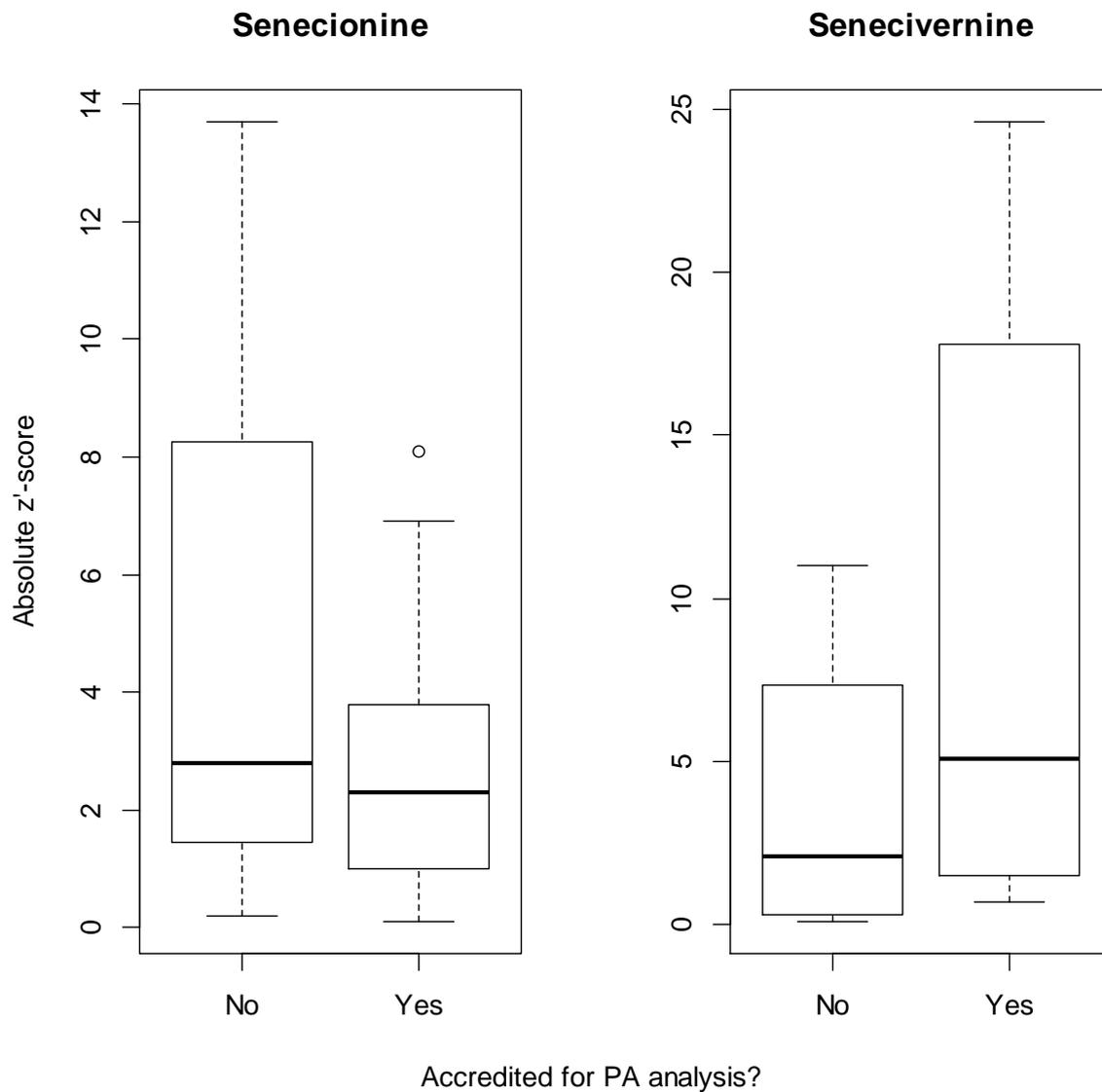


Figure 6- 11 Box & whisker plots of absolute z'-scores of Senecionine results in HO or Senecivernine results in HT grouped by whether an accreditation for PA analysis was held by the laboratories; bold horizontal line: median value, boxes extend from the first to the third quartile of the data.

Annex 7: Questionnaire data

Table 7- 1 Answers to questionnaire about experience and accreditation

| Lab Code | Previous experience | ISO 17025 accredited | Scope of accreditation | Lab Code | Previous experience | ISO 17025 accredited | Scope of accreditation |
|----------|---------------------|----------------------|--|----------|---------------------|----------------------|---|
| 1 | >12 months | Yes | DAkKS | 14 | >12 months | No | |
| 2 | >12 months | No | | 15 | >12 months | Yes | Honey, herbs and tea, grains |
| 3 | >12 months | No | | 16 | >12 months | No | |
| 4 | >12 months | Yes | tea, honey | 17 | 3-12 months | No | |
| 5 | 0-3 months | No | | 18 | 0-3 months | No | |
| 6 | >12 months | Yes | flexible scope, determination of contaminants, plant toxins and mycotoxins by LC-MS/MS | 19 | >12 months | No | |
| 7 | >12 months | Yes | Food | 20 | >12 months | Yes | DAkKS Number D-PL-XXXXX-01-00, flexible |
| 8 | >12 months | No | | 21 | >12 months | Yes | See on the DAKKS website under our accreditation number: D-PL-XXXXX-01-00 |
| 9 | >12 months | Yes | HPLC-MS/MS | 22 | >12 months | Yes | Bee Products, Honey, Plant materials, tea and herbal tea, spices, feed |
| 10 | >12 months | Yes | food (tea, herbal tea, spices) | 23 | >12 months | No | |
| 11 | >12 months | Yes | | 24 | 3-12 months | Yes | D-PL-XXXXX-01-00 |
| 12 | >12 months | Yes | Tea, herbal Tea, Honey | 25 | >12 months | Yes | HONEY |
| 13 | >12 months | Yes | flexible accreditation | 26 | >12 months | No | |

Table 7- 2 Answers to questionnaire about treatment of test item HO (Labs 1-13)

| Lab Code | HO test portion size [g] | HO extractant | HO extractant simplified | HO extractant vol. [mL] | HO extraction time [min] | HO cleaned up? | HO Clean-up type | HO calibration mode |
|----------|--------------------------|--|--------------------------|-------------------------|--------------------------|----------------|--|-----------------------------|
| 1 | 5 | water/acetonitrile | aqu/org | 20 | 30 | Yes | QuEChERS MIX I: MgSO ₄ , NaCl, Na ₂ H Citrate, Na ₃ Citrate; QuEChERS MIX III: Diamino, MgSO ₄ | |
| 2 | 0 | N/A | | | | No | | N/A |
| 3 | 0.2 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 5 | 720 | Yes | SCX-SPE | isotope dilution |
| 4 | 2 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 10 | 15 | Yes | SPE | matrix-matched external STD |
| 5 | 1.5 | H ₂ SO ₄ 0,1M and ACN | aqu/org/acidic | 10 | 30 | Yes | QUECHERS | neat solvent external STD |
| 6 | 5 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 30 | 10 | Yes | HR-X (Macherey Nagel, C18) | matrix-matched external STD |
| 7 | 10 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 30 | 30 | Yes | SPE with HF Bond Elut LRC-SCX | matrix-matched external STD |
| 8 | 10 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 30 | 30 | Yes | SPE cartridge SCX 500 mg | matrix-matched external STD |
| 9 | 2 | Water, H ₂ SO ₄ | aqu/acidic | 40 | 20 | Yes | filter | neat solvent external STD |
| 10 | 5 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 50 | 15 | Yes | SPE with Chromabond HR-X | standard addition |
| 11 | 2 | 2 % formic acid | aqu/acidic | 40 | 30 | Yes | centrifugated, filtered, SPE clean-up | matrix-matched external STD |
| 12 | 10 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 30 | 30 | Yes | Cation exchange (HF Bond Elut LRC-SCX, 500 mg) | matrix-matched external STD |
| 13 | 5 | acetonitrile (adjusted to alkaline pH w/ ammonia solution) | aqu/org/basic | 10 | 10 | No | | matrix-matched external STD |

Table 7- 3 Answers to questionnaire about treatment of test item HO (Labs 14-26)

| Lab Code | HO test portion size [g] | HO extractant | HO extractant simplified | HO extractant vol. [mL] | HO extraction time [min] | HO cleaned up? | HO Clean-up type | HO calibration mode |
|----------|--------------------------|---|--------------------------|-------------------------|--------------------------|----------------|---|---|
| 14 | 5 | acetonitrile:water, 1:1 (QuEChERS-like method) | aqu/org | 20 | 2 | No | | matrix-matched external STD |
| 15 | 0.5 | Acetonitrile : Methanol : Water 1:1:1 (v:v:v) | aqu/org | 10 | 15 | No | | neat solvent external STD / standard addition |
| 16 | 5 | 20 g NaCl in 800 ml H ₂ O, making up to 1 L with HCl 37% | aqu/acidic | 20 | 60 | No | | neat solvent external STD |
| 17 | 10 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 30 | 30 | Yes | MCX SPE_Waters | Prespiked matrix matched calibration |
| 18 | 1 | 0.1 n HCl in Water | aqu/acidic | 20 | 60 | Yes | NH ₂ | matrix-matched external STD |
| 19 | 4 | Hydrochloric acid and Acetonitrile | aqu/org/acidic | 8 | 10 | Yes | MgSO ₄ and PSA | matrix-matched external STD |
| 20 | 1 | aqueous citric buffer solution | aqu/acidic | 25 | 15 | Yes | Solid phase extraction | matrix-matched external STD |
| 21 | 10 | no further information | | 0 | 0 | Yes | no further information | standard addition |
| 22 | 1 | Acetonitrile | org | 4 | 6 | Yes | liquid liquid extraction | standard addition |
| 23 | 10 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 30 | 30 | Yes | SPE Bond Elut LRC SCX 500mg | matrix-matched external STD |
| 24 | 0 | not | | 0 | 0 | No | | |
| 25 | 20 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 200 | 30 | Yes | SCX SPE CARTRIDGE | standard addition |
| 26 | 2 | 0.2% formic acid in water | aqu/acidic | 20 | 30 | Yes | Solid Phase Extraction on Phenomenex Strata X 200 mg/6 ml | standard addition |

Table 7- 4 Answers to questionnaire about treatment of test item HT (Labs 1-13)

| Lab Code | HT test portion size [g] | HT extractant | HT extractant simplified | HT extractant vol. [mL] | HT extraction time [min] | HT cleaned up? | HT Clean-up type | HT calibration mode |
|----------|--------------------------|--------------------------|--------------------------|-------------------------|--------------------------|----------------|---|-----------------------------|
| 1 | 0 | not analyzed | | 0 | 0 | No | | N/A |
| 2 | 2 | acetic acid and methanol | aqu/org/acidic | 20 | 23 | No | | standard addition |
| 3 | 0.1 | 0.05 mol/l H2SO4 | aqu/acidic | 5 | 720 | Yes | SCX-SPE | isotope dilution |
| 4 | 10 | 0.05 mol/l H2SO4 | aqu/acidic | 200 | 15 | Yes | SPE | standard addition |
| 5 | 1 | 0.05 mol/l H2SO4 | aqu/acidic | 25 | 30 | Yes | SPE MCX | neat solvent external STD |
| 6 | 2 | 0.05 mol/l H2SO4 | aqu/acidic | 40 | 30 | Yes | HR-X (Macherey Nagel, C18) | standard addition |
| 7 | 2 | 0.05 mol/l H2SO4 | aqu/acidic | 40 | 30 | Yes | SPE with DSC-C18 SPE (Supelco) | neat solvent external STD |
| 8 | 2 | 0.05 mol/l H2SO4 | aqu/acidic | 40 | 30 | Yes | SPE DSC C18 | standard addition |
| 9 | 2 | Water, H2SO4 | aqu/acidic | 40 | 20 | Yes | filter | neat solvent external STD |
| 10 | 2.5 | 0.05 mol/l H2SO4 | aqu/acidic | 50 | 30 | Yes | SPE with Chromabond HR-X | standard addition |
| 11 | 2 | 2 % formic acid | aqu/acidic | 40 | 30 | Yes | centrifugated, filtered, SPE clean-up | matrix-matched external STD |
| 12 | 2 | 0.05 mol/l H2SO4 | aqu/acidic | 20 | 15 | Yes | C18-SPE (Discovery DSC18, 6 ml, 500 mg) | neat solvent external STD |
| 13 | 2 | 0.05 mol/l H2SO4 | aqu/acidic | 40 | 15 | Yes | SPE (C18-ec) | matrix-matched external STD |

Table 7- 5 Answers to questionnaire about treatment of test item HT (Labs 14-26)

| Lab Code | HT test portion size [g] | HT extractant | HT extractant simplified | HT extractant vol. [mL] | HT extraction time [min] | HT cleaned up? | HT Clean-up type | HT calibration mode |
|----------|--------------------------|---|--------------------------|-------------------------|--------------------------|----------------|--|--|
| 14 | 2 | acetonitrile:water, 1:1 (QuEChERS-like method) | aqu/org | 20 | 30 | No | | matrix-matched external STD |
| 15 | 0.5 | Acetonitrile : Methanol : Water 1:1:1 (v:v:v) | aqu/org | 10 | 15 | No | | neat solvent external STD / standard addition |
| 16 | 2.5 | 20 g NaCl in 800 ml H ₂ O, making up to 1 L with HCl 37% | aqu/acidic | 25 | 60 | No | | neat solvent external STD |
| 17 | 2 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 50 | 30 | Yes | MCX SPE waters | Prespiked matrix matched calibration |
| 18 | 0.5 | 0.1 n HCl in Water | aqu/acidic | 20 | 60 | Yes | GCB | matrix-matched external STD |
| 19 | 4 | Hydrochloric acid and Acetonitrile | aqu/org/acidic | 8 | 10 | No | | NA |
| 20 | 1 | aqueous citric buffer solution | aqu/acidic | 25 | 15 | Yes | Solid phase extraction | matrix-matched external STD |
| 21 | 1 | no further information | | 0 | 0 | Yes | no further information | standard addition |
| 22 | 2 | no further information | | 0 | 0 | Yes | no further information | standard addition |
| 23 | 2 | 0.05 mol/l H ₂ SO ₄ | aqu/acidic | 40 | 30 | Yes | SPE DSC-C18 (Supelco) 500mg | matrix-matched external STD |
| 24 | 2 | Water+sulfuric acid | aqu/acidic | 20 | 15 | Yes | SPE | |
| 25 | 1 | 0.2% formic acid in water | aqu/acidic | 40 | 30 | Yes | Solid Phase Extraction on Phenomenex Strata X 200 mg/6 ml | standard addition |
| 26 | 2 | 0.2% formic acid in water | aqu/acidic | 40 | 30 | Yes | Solid Phase Extraction on Phenomenex Strata X 200 mg/6 ml | standard addition |

Table 7- 6 Answers to questionnaire about separation conditions (Labs 1-13)

| Lab Code | Separation method | Sep. meth. other | Analytical column | Mobile phase | Flow rate [ml/min] |
|----------|-----------------------|------------------|---|---|--------------------|
| 1 | Liquid chromatography | | C18, 50 mm, 2,0 mm, 2µm (YMC Ultra Hydrosphere 50*2,0mm; (Partikel 2µm) (Poren 12mm) | eluent A: distilled water (1 %CH ₃ COOH, 0,05 % 10 M NH ₄ (HCOO)) ; eluent B: methanol (1 %CH ₃ COOH, 0,05 % 10 M NH ₄ (HCOO)) | 0.2 |
| 2 | Liquid chromatography | | Thermo Hypersil Gold C18; 150x2,1mm | water and methanol | 0.3 |
| 3 | Liquid chromatography | | Synergi Max-RP 4µm, 150 x 2 mm | H ₂ O/ACN Gradient, each phase containing 0.3 % formic acid | 0.3 |
| 4 | Liquid chromatography | | C18, 150 x 2.0 mm | water/acetonitrile with ammonium formiate and formic acid | 0.3 |
| 5 | Liquid chromatography | | Luna C8 150 x 2 mm 3 µm | ammonio formiate in H ₂ O and ammonio formiate in MeOH | 0.25 |
| 6 | Liquid chromatography | | C18, 100mm, 2,1mm, 1,8 µm (Eclipse Plus) | A: Methanol+0,1% FA B: Wasser+0,1% FA | 0.3 |
| 7 | Liquid chromatography | | Kinetex C18, 100 mm, 2.1 mm, 2,6 µm | 5 mM NH ₄ COO and 0.1% COOH in water/ 5 mM NH ₄ COO and 0.1% COOH in methanol | 0.2 |
| 8 | Liquid chromatography | | C18, 200 mm, 2.1 mm, 1.9 µm, Hypersil Gold | water / Methanol with HCOOH 0.1% and ammonium formiate 315 mg/l | 0.25 |
| 9 | Liquid chromatography | | N/A | N/A | 0.4 |
| 10 | Liquid chromatography | | Acquity UPLC BEH C18 1,7 µm, 2,1 x 100 mm | ammonium formiate + formic acid in water / methanol | 0.2 |
| 11 | Liquid chromatography | | Waters, ACQUITY UPLC BEH C18 150 x 2,1 mm 1,7 µm | 0,1 % formic acid in water; 0,1 % formic acid in acetonitrile | 0.2 |
| 12 | Liquid chromatography | | RP18 (Kinetex 2.6 µ EVO C18 100A, 150 * 2.1 mm) | Gradient H ₂ O and MeOH (both 5 mM NH ₄ HCO ₂ and 0.1% CH ₂ O ₂) | 0.35 |
| 13 | Liquid chromatography | | Waters BEH C18 (1.7 µm, 100 mm, 2.1 mm ID) | A: 0.005 M ammonium formate / formic acid in water (pH 3); B: A: 0.005 M ammonium formate / formic acid in methanol (pH 3) | 0.35 |

Table 7- 7 Answers to questionnaire about separation conditions (Labs 14-26)

| Lab Code | Separation method | Sep. meth. other | Analytical column | Mobile phase | Flow rate [ml/min] |
|----------|-----------------------|------------------------|---|---|--------------------|
| 14 | Liquid chromatography | | silica-based reversed-phase (C18) analytical column Acquity UPLC HSS T3 (length 100 mm, 2.1 mm inner diameter, 1.8 µm particle size) | 5mM ammonium formate and 0.2% formic acid both in methanol (A) and milli-Q water (B) | 0.3 |
| 15 | Liquid chromatography | | Waters Acquity UPLC BEH C18 1.7µm 2.1 x 150 mm | water : methanol + ammoniumformiat and formic acid | 0.3 |
| 16 | Liquid chromatography | | Phenomenex Kinetex 2.6 µm XB-C18 column (100 x 4.6 mm) | MPA: H2O + 5 mM ammonium acetate + 0.05% acetic acid; MPB: acetonitrile | 1.5 |
| 17 | Other | LCMS-MS | C18,100mm,2.1um,3um | 5mM ammonium acetate ,0.1% formic acid in MeOH/water | 0.4 |
| 18 | Liquid chromatography | | Agilent InfinityLab Poroshell 120 EC-C18,2.1 x 100 mm, 2.7 µm (p/n 695775-902) with HPLC, Guard (p/n 821725-911) at 25 °C | 0.025 % formic acid in 5 mM ammonium formate in Water (A), 0.025 % formic acid 5 mM ammonium formate in methanol (B), | 0.3 |
| 19 | Liquid chromatography | | 150mm*2mm, C18 | Methanol/ammonium acetate | 0.15 |
| 20 | Liquid chromatography | | Thermo Scientific TM Hypersil Gold C18 150 x 3 mm, 3 µm particle size | A: methanol + 5 mmol ammonium formate + 1% formic acid; B:water + 5 mmol ammonium formate + 1% formic acid | 0.4 |
| 21 | Liquid chromatography | no further information | no further information | no further information | 0 |
| 22 | Liquid chromatography | | C18, 150mm, 2.1mm, 5µm | formic acid in water / formic acid in methanol | 0.3 |
| 23 | Liquid chromatography | | Acclaim Vanquish PA2, 150 x 2.1 mm, 2.2 µm particle size | A: Ammonium formate (315mg/L): formic acid (1ml/L) in water; B: Ammonium formate (315mg/L): formic acid (1ml/L) in methanol | 0.3 |
| 24 | Other | LC MS/MS | Hypersil gold 150 x 2,1mm; 175A, 3µm | Eluent A: Water Eluent B: Methanol | 0.3 |
| 25 | Liquid chromatography | | C18, 150MM, 2.1MM,1.9UM | A:5mM AMMONIUM FORMATE, 0.1% FORMIC ACID IN METHANOL; B: 5mM AMMONIUM FORMATE, 0.1% FORMIC ACID IN ULTRAURE WATER | 0.3 |

Table 7- 8 Answers to questionnaire about detection conditions (Labs 1-13)

| Lab Code | Detector | Ionization mode | Ionization mode other | MS mode | PAs determined as | PANOs determined as | Other | Derivatization? | Which deriv. |
|----------|-------------------|-----------------|-----------------------|-------------|---|--|-------|-----------------|------------------------|
| 1 | Mass spectrometry | ESI | | MRM | Individual PAs | measured as such | | No | |
| 2 | Mass spectrometry | ESI | | SRM | Individual PAs | measured as such | | No | |
| 3 | Mass spectrometry | ESI | | MRM | Sum after conversion of PAs to common analyte | reduced to the amine before separation | | Yes | Phthalic acid anhydrid |
| 4 | Mass spectrometry | ESI | | SRM | Individual PAs | measured as such | | No | |
| 5 | Mass spectrometry | ESI | | SRM | Individual PAs | measured as such | | No | |
| 6 | Mass spectrometry | ESI | | MRM | Individual PAs | measured as such | | No | |
| 7 | Mass spectrometry | ESI | | sMRM | Individual PAs | measured as such | | No | |
| 8 | Mass spectrometry | ESI | | MRM | Individual PAs | measured as such | | No | |
| 9 | Mass spectrometry | ESI | | | Individual PAs | measured as such | | No | |
| 10 | Mass spectrometry | ESI | | MRM | Individual PAs | measured as such | | No | |
| 11 | Mass spectrometry | ESI | | MRM | Individual PAs | measured as such | | No | |
| 12 | Mass spectrometry | ESI | | MRM | Individual PAs | measured as such | | No | |
| 13 | Mass spectrometry | ESI | | dynamic MRM | Individual PAs | measured as such | | No | |

Table 7- 9 Answers to questionnaire about detection conditions (Labs 14-26)

| Lab Code | Detector | Ionization mode | Ionization mode other | MS mode | PAs determined as | PANOs determined as | Other | Derivatization? | Which deriv. |
|----------|-------------------|-----------------|------------------------|------------------------|-------------------|---------------------|---|-----------------|--------------|
| 14 | Mass spectrometry | ESI | n | HRMS/MS | Individual PAs | measured as such | | No | |
| 15 | Mass spectrometry | ESI | n | MRM | Individual PAs | measured as such | | No | |
| 16 | Mass spectrometry | ESI | n | SRM | Individual PAs | measured as such | | No | |
| 17 | Mass spectrometry | ESI | n | MRM | Individual PAs | measured as such | | No | |
| 18 | Mass spectrometry | ESI | n | SRM | Individual PAs | measured as such | | No | |
| 19 | Mass spectrometry | Other | EI | SIM | Other | measured as such | some as individual, some as sum of several PA (eg. Lycopsamin/Indicin/Intermedin) | No | |
| 20 | Mass spectrometry | ESI | n | SIM | Individual PAs | measured as such | | No | |
| 21 | Mass spectrometry | Other | no further information | no further information | Individual PAs | measured as such | | TRUE | |
| 22 | Mass spectrometry | ESI | n | SRM | Individual PAs | measured as such | | No | |
| 23 | Mass spectrometry | ESI | ESI positive mode | SIM | Individual PAs | measured as such | | No | |
| 24 | Mass spectrometry | ESI | n | SRN | Individual PAs | measured as such | | No | |
| 25 | Mass spectrometry | ESI | n | MRM | Individual PAs | measured as such | | No | |
| 26 | Mass spectrometry | ESI | n | MRM | Individual PAs | measured as such | | No | |

Table 7- 10 Comments from Labs 1-13

| Lab Code | Comments |
|----------|---|
| 1 | |
| 2 | |
| 3 | |
| 4 | |
| 5 | I haven't the Indicine and Integerrimine standards, so I could't resolve these PAs from its isomers (Intermedine/Lycopsamine and Senecionine/Senecivernine). The sum of all PAs is the sum of the only PAs with a concentration >LOQ (the other PAs are considered as 0 ppb) |
| 6 | Indicine and Integerrimine are not part of our method |
| 7 | The Concentration and composition of the honey sample did not match our common samples (Concentration too high and too much N-Oxides) |
| 8 | indicine or intermedine are not separated on a C18 column |
| 9 | N/A |
| 10 | |
| 11 | On mass transiotion of Senecionone-NO we recognized a Peak with shifted retention time to Senecionine NO Peak. It also has differing ion rations. The Peak is poparbly a different PA for which we do not ahve any standard. The transiotions are typical for PA. Results of Senecivernine adn Senecionine are reported as the sum of both. |
| 12 | N/A |
| 13 | Wouldn't it have been wiser to report individual PA results and summarize them subsequently? |

Table 7- 11 Comments from Labs 14-26

| Lab Code | Comments |
|----------|--|
| 14 | Under conditions described within the PT, several PAs were difficult to resolve chromatographically and spectrometrically (lycopsamine x indicine x intermedine; senecionine x integerrimine x senecivernine; retrorsine x jacobine); the reported results might be therefore affected by the co-occurrence of these PAs. |
| 15 | Low amounts of Senecivernin and the respective N-oxide, just below the respective LOQs where found in the herbal tea. Reporting as the sum of amine and the respective N-oxide is disadvantageous to accuracy and comparability. |
| 16 | The results are reported as sum of PA and its corresponding N-oxide, however, echimidine-N-oxide, intermedine-N-oxide, lycopsamine-N-oxide were not measured. Under our LC conditions there is co-elution of lycopsamine and indicine. HT material: a peak interfering with senecionine-N-oxide prevented proper peak integration, therefore result for this PANO is not reported. |
| 17 | |
| 18 | LOQs: very strong dependence on matrix components |
| 19 | |
| 20 | Meaning of "coverage factor" is unclear. |
| 21 | Section 3.3., 3.4., 4.3., 4.4., 5.4., 5.5.2 no further information |
| 22 | |
| 23 | Methods used were the BfR methods for honey and plant material. Tea sample was extracted as 2 x 20ml extraction solvent for 15 min each extraction. |
| 24 | only tea, sum of amine and N-oxide |
| 25 | ION RATIO DISTORTION OBSERVED FOR RETRORSINE-N-OXIDE IN TEA SO VALUE FOR RETRORSINE NOT REPORTED. ION RATION DISTORTION WAS ALSO OBSERVED FOR ECHIMIDINE IN HONEY. |
| 26 | Quantification was performed by means of standard addition to the sample and using internal standard correction |

Annex 8: Participating laboratories

| Organisation | Part. key | Firstname | Surname | Country |
|---|-------------|------------|--------------|----------------|
| Agri-Food and Veterinary Authority of Singapore | PLA11131621 | Lim | Poh Choo | SINGAPORE |
| Austrian Agency for Health and Food Safety | CCAK1172028 | Christoph | Czerwenka | AUSTRIA |
| Bayerisches Landesamt für Gesundheit und Lebensmittelsicherheit | GFBS1348628 | Franziska | Gaßmann | GERMANY |
| Chemical and Veterinary Analytical Institute Muensterland-Emscher-Lippe | KOC41110468 | Oliver | Keuth | GERMANY |
| Chemisches und Veterinäruntersuchungsamt Ostwestfalen-Lippe | AEC41082382 | Elisabeth | Apel | GERMANY |
| CVUA Stuttgart | KTC10323649 | Thomas | Kapp | GERMANY |
| Eurofins Dr. Specht Laboratorien GmbH | SEEQ1345165 | Edda | Sassen | GERMANY |
| Fera Science Ltd. | MSFF1153945 | Susan | MacDonald | UNITED KINGDOM |
| FoodQS GmbH | RNF12432073 | Nadine | Raum | GERMANY |
| Gesellschaft für Bioanalytik mbH | JSG11185486 | Stefan | Jäger | GERMANY |
| HEALTH SCIENCES AUTHORITY | LAHF1128167 | Angela | Li | SINGAPORE |
| Intertek Food Services GmbH Bremen | ZHI12913009 | Hauke | Zinow | GERMANY |
| IZSLER | CEIC1078150 | ELISABETTA | CAPRAI | ITALY |
| Landesbetrieb Hessisches Landeslabor (LHL) | KJL13694049 | Johannes | Kemme | GERMANY |
| Landesuntersuchungsamt Rheinland-Pfalz | KSLI1049678 | Simone | Kasper | GERMANY |
| Lower Saxony State Office for Consumer Protection and Food Safety | DKL12243546 | Kay | Dietrichkeit | GERMANY |
| PiCA GmbH | RAPF1342087 | Anna | Romanotto | GERMANY |
| PMA Sindelfingen GmbH | GSP10165902 | Stefan | Glöckler | GERMANY |
| QSI - Quality Services International GmbH & Co.KG | WTQQ1279373 | Tobias | Wiezorek | GERMANY |
| RIKILT Wageningen University & Research | MHRN1018129 | Hans | Mol | NETHERLANDS |
| SCL laboratoire de Strasbourg | GLSV1151252 | Ledoux | Gérald | FRANCE |
| SGS Germany GmbH | MNS12889924 | Nicolaus | Mouillard | GERMANY |
| Teekanne GmbH & Co.KG | BATL1197037 | Andreas | Bollert | GERMANY |
| TU Braunschweig | BTTI1231664 | Till | Beuerle | GERMANY |
| University of Chemistry and Technology Prague | KVUF1247054 | Vladimir | Kocourek | CZECH REPUBLIC |
| Veterinary and Agrochemical Research Centre (CODA-CERVA) | MSVT1247439 | Svetlana | Malysheva | BELGIUM |

GETTING IN TOUCH WITH THE EU

In person

All over the European Union there are hundreds of Europe Direct information centres. You can find the address of the centre nearest you at: <http://europa.eu/contact>

On the phone or by email

Europe Direct is a service that answers your questions about the European Union. You can contact this service:

- by freephone: 00 800 6 7 8 9 10 11 (certain operators may charge for these calls),
- at the following standard number: +32 22999696, or
- by electronic mail via: <http://europa.eu/contact>

FINDING INFORMATION ABOUT THE EU

Online

Information about the European Union in all the official languages of the EU is available on the Europa website at: <http://europa.eu>

EU publications

You can download or order free and priced EU publications from EU Bookshop at: <http://bookshop.europa.eu>. Multiple copies of free publications may be obtained by contacting Europe Direct or your local information centre (see <http://europa.eu/contact>).

JRC Mission

As the science and knowledge service of the European Commission, the Joint Research Centre's mission is to support EU policies with independent evidence throughout the whole policy cycle.



EU Science Hub
ec.europa.eu/jrc



@EU_ScienceHub



EU Science Hub - Joint Research Centre



Joint Research Centre



EU Science Hub

